
Gebruiksaanwijzing - Draadloze Infrarood Rugband
User Manual - Wireless Infrared Back Belt
Manual d'utilisation - Ceinture dorsale infrarouge sans fil
Gebrauchsanweisung - Kabelloser Infrarot-Rückengurt

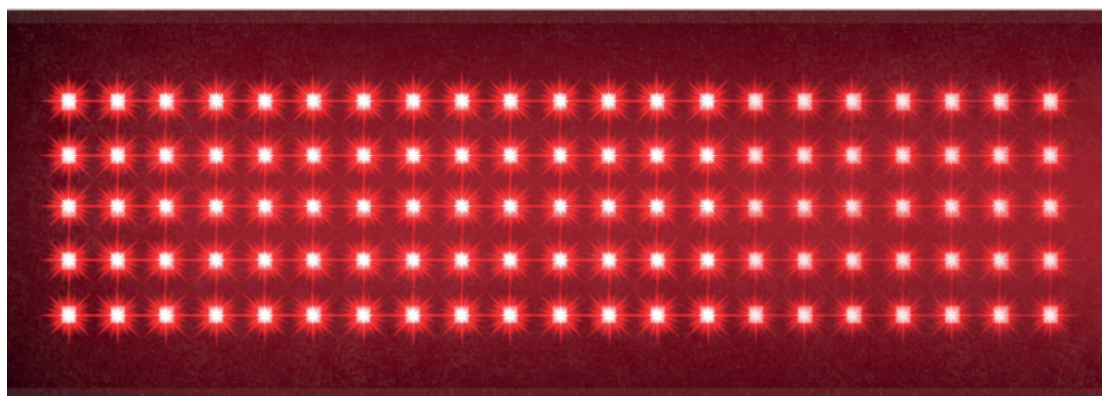Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using 'the Wireless Infrared Back Belt ', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it!
Read the instructions for use carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- This product is not suitable for individuals (including children) with physical or mental limitations, or for those who lack experience or knowledge of this product. Use of the product is only permitted if supervised by a responsible person who ensures the user's safety. Keep and store this product out of the reach of children. We recommend using this product only if you are 16 years of age or older.
- Do not use this product in combination with implanted electronic devices in the body, such as pacemakers, artificial joints, etc. Consult your doctor if in doubt.
- Do not use this product if you are pregnant, have a severe heart condition, are seriously ill, or are an emergency patient. Consult your doctor if unsure.
- Do not use this product near babies.
- Do not use this product while sleeping.
- Do not use this device for more than 30 minutes in one therapy session.
- Do not use this product on open or fresh wounds.
- If you lend this product to someone else, make sure to provide the instruction manual with it.

Elektra

- Do not disassemble, repair, or modify the product. Doing so may cause damage. Warranty is void if you open or repair the product yourself. For safety, parts may only be replaced by authorized personnel. Always return the product to a service center for repairs to prevent damage or personal injury.
- Do not use the product near a bath, shower, swimming pool, or sink filled with water.
- Ensure the product does not come into contact with water. Do not immerse in water. If water does enter the product, switch it off immediately and contact Customer Service.
- Do not touch the product with wet hands.

Other

- Ensure that the buttons are not pressed unintentionally when plugging in or moving the product.
- Do not insert objects into the openings or slots of the product. Ensure that no foreign objects enter the product housing.

About the Wireless Infrared Back Belt

The Wireless Infrared Back Belt uses a combination of red light (660 nm) and near-infrared light (850 nm). These specific wavelengths are considered ideal for supporting recovery, skin improvement, and overall health.

Light therapy lamps are designed to stimulate the mitochondria – the energy-producing units of our cells. By improving mitochondrial function, the overall health of the body can be enhanced. Red light reaches depths that UV, blue, and green light cannot, and has a healing effect on the skin. Within the infrared spectrum, the rule applies: the shorter the wavelength, the deeper the radiation penetrates. Near-infrared light (short wave) is barely visible but penetrates the deepest, reaching muscles, joints, and deeper tissues more effectively.

Recommended Use

- Before use, check that your blood pressure is normal and do not wear accessories such as belts or jewelry under the device.
- Drink 200–300 ml of water (optionally with a little salt) before, during, and after each session to keep the body hydrated and flush out toxins.
- For first-time use, it is recommended to select the low or medium temperature setting.
- If possible, use the device after taking a bath for better effect.
- Do not use the device within one hour before or after meals. Avoid use when extremely hungry, with a full stomach, or when feeling very weak.
- After use, avoid exposure to air conditioning and do not bathe within one hour.
- Light therapy works best when applied consistently. Best results are achieved with daily use, but it is recommended to use the device at least 4 times per week. For serious injuries, chronic pain, or severe skin conditions, using the device up to 3 times per day may be beneficial.
- Do not use the device for more than 30 minutes at a time and wait at least 2 hours between 30-minute sessions.
- **Always consult a doctor before use in case of contraindications, medical conditions, or when in doubt.**

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents

- Good4You Wireless Infrared Back Belt- Model FE-7509

Included Accessories

- 1x User Manual
- 1x Remote Control
- 1x USB Cable

Startup Instructions


Charging the device

- ① This product has a built-in rechargeable battery. For the first use, it is recommended to fully charge the battery.
1. Connect the USB cable to the Infrared Back Belt and to an external power source: via a USB adapter (5V 2A) in a wall socket or with a power bank. During charging, the battery icon will slowly flash with a white light.
 2. Once the battery is fully charged, the battery indicator on the remote control will light up white. The estimated charging time is 2.5 hours.
 3. When the battery is almost empty, the battery indicator will turn red and the device will shut down automatically. Recharge the device before using it again.

Using the device

1. Wear the Infrared Back Belt around your waist and place it in a comfortable position. (Note: the belt can also be used on your back, shoulders, arms, legs, hips, or abdomen.)
2. Press the "  " button for 1.5 seconds to turn the Infrared Back Belt on or off. The default setting is: **Mode 1 / L3 / Timer 15** minutes.
3. Short press the "  " button to set the usage time. There are 4 time settings: **15, 20, 25, or 30 minutes.**
4. Press the "  " button to adjust the light intensity. There are 4 intensity levels: **L1, L2, L3, L4.**
5. Press the " **M** " button to set the mode. There are 4 modes:
 - **Mode 1:** 850Hz + 660Hz
 - **Mode 2:** 660Hz
 - **Mode 3:** 850Hz
 - **Mode 4:** 10Hz
6. Press the "  " button to set the vibration. There are 4 vibration levels: **M1, M2, M3, M4.**

Remote Control

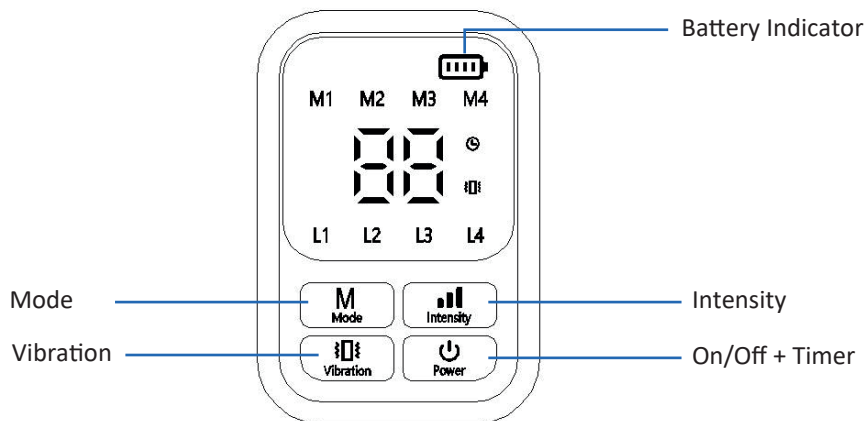Technical data

Power	12W
Battery Capacity	5000 mAh
Number of LEDs	105 pcs
Product dimensions	130 x 18 x 11 cm

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.

- Turn off and unplug the product before cleaning. Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Extend warranty for FREE

1) Scan the QR code below to find Good4You's Facebook page and then 'Like' it.

2) Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 1 year to 3 years.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

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© OnlineGoods