
Gebruiksaanwijzing - Draagbaar Roodlichttherapie-apparaat
User Manual - Portable Red Light Therapy Device
Manual d'utilisation - Appareil portable de lumière rouge
Gebrauchsanweisung - Tragbares Rotlichttherapiegerät

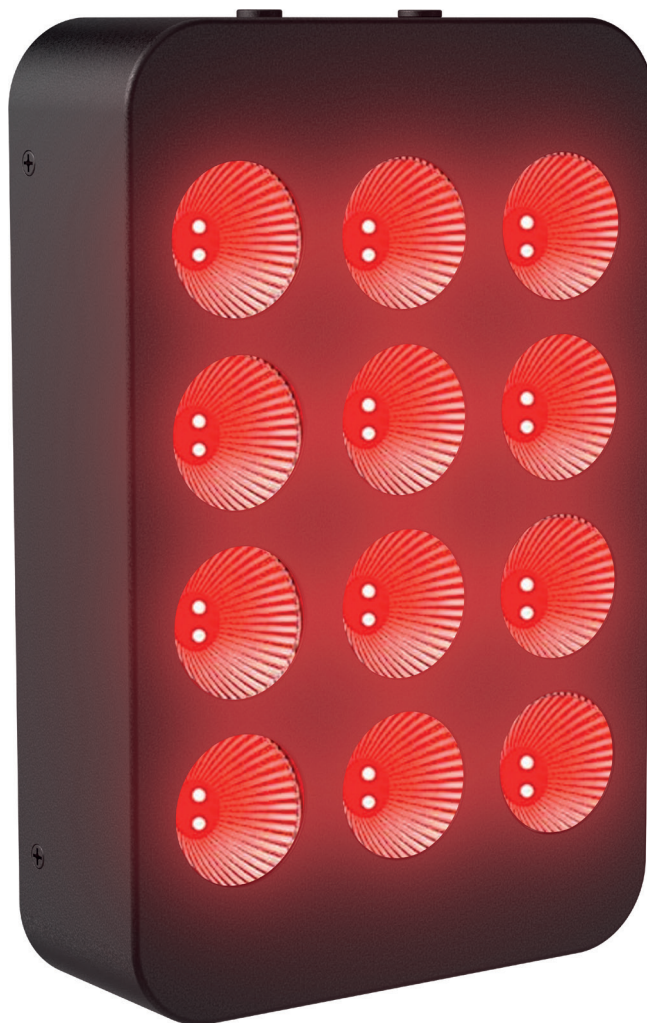Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using 'the Portable Red Light Therapy Device', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it!
Read the instructions for use carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- This product is not suitable for individuals (including children) with physical or mental limitations, or for those who lack experience or knowledge of this product. Use of the product is only permitted if supervised by a responsible person who ensures the user's safety. Keep and store this product out of the reach of children. We recommend using this product only if you are 16 years of age or older.
- Do not use this product in combination with implanted electronic devices in the body, such as pacemakers, artificial joints, etc. Consult your doctor if in doubt.
- Do not use this product if you are pregnant, have a severe heart condition, are seriously ill, or are an emergency patient. Consult your doctor if unsure.
- Do not use this product near babies.
- Do not use this product while sleeping.
- Do not use this device for more than 20 minutes in one therapy session.
- Medicines, cosmetics, or dietary supplements may cause hypersensitivity or allergic reactions on the skin. Stop the treatment immediately if this occurs.
- Do not use this product on open or fresh wounds.
- Do not use this product in combination with ointments, balms, or creams that produce heat. This may cause skin burns.
- If you lend this product to someone else, make sure to provide the instruction manual with it.

Elektra

- Do not disassemble, repair, or modify this product yourself. Doing so may damage the product. Warranty will be void if the product is opened or repaired by unauthorized persons. For your safety, such repairs must only be performed by authorized service staff. Always return the product to customer service to avoid damage or personal injury.
- Do not use this product near a bath, shower, swimming pool, or any other container filled with water.
- Ensure this product does not come into contact with water. Do not immerse the product in water. If water does get into the product, stop using it immediately and contact customer service.
- Do not touch the product with wet hands.
- Do not cover the product while it is in use. Never place pillows, blankets, or other objects on top of the device as this may cause fire, electric shock, personal injury, or damage.
- Ensure ventilation openings remain unobstructed during use.

Other

- Ensure that the buttons are not pressed unintentionally when plugging in or moving the product.
- Do not insert objects into the openings or slots of the product. Ensure that no foreign objects enter the product housing.

Recommended Use

- The ideal usage frequency is between 3 and 7 times per week, maximum once per day per body part.
- Start slowly: position the lamp only for a few minutes at a greater distance so your body can get used to it.
- **Guideline:** the closer the lamp is to the skin, the shorter the recommended treatment duration. Recommended duration is 5 to 20 minutes at a distance of up to 90 cm.
- Example of a recommended treatment: 12 minutes at 30 cm distance.
- Use the lamp for a maximum of 20 minutes. If you want to use the lamp for 20 minutes, increase the distance.
- The lamp has 3 modes: RED infrared light (660 nm), NIR infrared light (850 nm), or a combination of both.
- Note! The lamp **does not** emit heat. The lamp is designed exclusively for LED-light therapy.
- Start the desired program, then put on your protective glasses and try not to look directly into the lamp during treatment.
- Ensure the body part you want to treat is fully exposed to the lamp's light and not covered by clothing.
- **Always consult a doctor before use in case of contraindications, medical conditions, or if in doubt.**

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents


- Good4You Portable Red Light Therapy Device - Model AL60

Included Accessories

- 1x Travel case
- 1x User manual
- 1x USB cable

Control Panel

- | | |
|-----------------------|----------------------------|
| ①: Brightness control | ②: Digital display |
| ③: Charging light | ④: Battery indicator light |
| ⑤: ON/OFF | ⑥: USB interface |

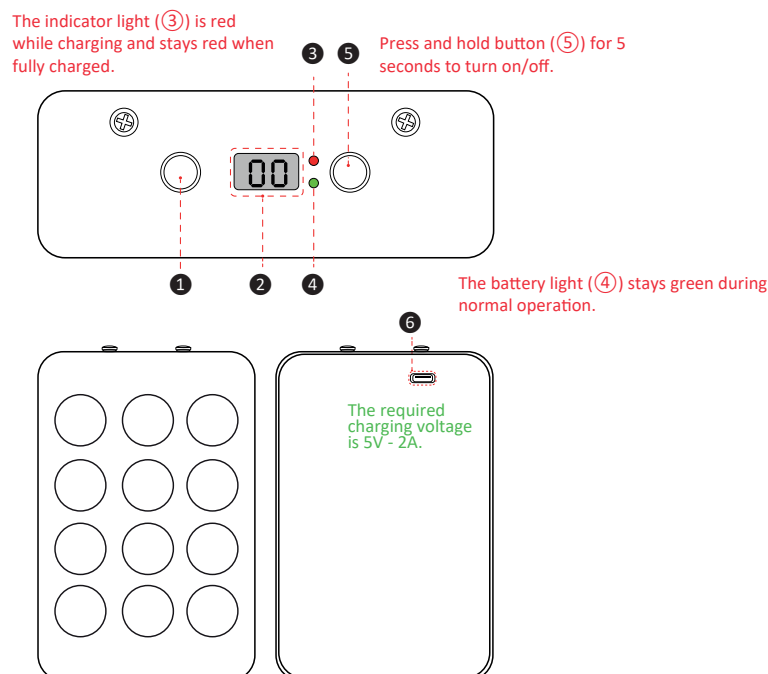Technical data

Power	60W
Battery	5000mA
Number of LEDs	12 (660nm 12 pcs + 850nm 12 pcs)
Product dimensions	13 × 8.5 × 3 cm
Product weight	480 g

Recommended use

The following guidelines are general advice for common health and wellness concerns. Your personal treatment needs may vary.

Skin health (wrinkles, acne, wounds)

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 3–5 times
- **Target area:** All skin areas with visible complaints
- **Distance from the Good4You lamp:** 10 cm

Joint and arthritis pain

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 6–7 times (if severe, possibly twice daily)
- **Target area:** Painful joint
- **Distance from the Good4You lamp:** 10 cm

Weight loss

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 5–7 times (best results before or after exercise)
- **Target area:** Area with cellulite or fat, or the entire body
- **Distance from the Good4You lamp:** 10 cm

Muscle recovery

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 3–5 times
- **Target area:** Stressed muscles (e.g., after sports)
- **Distance from the Good4You lamp:** 10 cm

Sleep and circadian rhythm

- **Treatment time per session:** 10–15 minutes
- **Sessions per week:** 5–7 times (preferably 30 minutes before sleeping)
- **Target area:** Whole body; place the lamp where you sleep
- **Distance from the Good4You lamp:** Free to determine in the bedroom

Brain health and cognitive support

- **Treatment time per session:** 10 minutes, twice daily
- **Sessions per week:** 5–7 times
- **Target area:** Head
- **Distance from the Good4You lamp:** 10 cm

Boosting Testosterone

- **Treatment time per session:** 5–10 minutes, 2× per day
- **Sessions per week:** 5–7 times
- **Target area:** Testicles
- **Distance from the Good4You lamp:** 10 cm

Thyroid Function

- **Treatment time per session:** 8–10 minutes
- **Sessions per week:** 5–7 times
- **Target area:** Thyroid area (slightly tilt head backward)
- **Distance from the Good4You lamp:** 10 cm

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.

- Turn off and unplug the product before cleaning. Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Extend warranty for FREE

1) Scan the QR code below to find Good4You's Facebook page and then 'Like' it.

2) Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 1 year to 3 years.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

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© OnlineGoods