Gebruiksaanwijzing - Infrarood Rugband
User Manual - Infrared Backband
Manual d'utilisation - Bande arrière infrarouge
Gebrauchsanweisung - Infrarot-Backband

Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using 'the infrared belt', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it! Read the instructions for use carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or people who have no experience with or knowledge of the product. Use of the product is permitted only when a supervisor is present and responsible for the safety of the user. Use and store the product out of the reach of children. We recommend using the product only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; are pregnant; have injuries or lesions to your foot(s), or have conditions such as warts, swimmer's eczema, rashes, or open wounds; have a malignant tumor, heart disease, osteoporosis, or a broken spine; have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nervous disorders, acute inflammation, skin diseases, broken skin, or pain without obvious or immediate cause.
- Do not use the product with more than one person at a time. Do not use the product near your eyes or other sensitive areas. Do not use the product when you are sleeping. Do not use the product when you are driving a vehicle or operating machinery.
- If you lend the product to someone else, be sure to give the product along with the manual.
- Prolonged use of the product may result in overheating. After 10 minutes of use, take a break for at least 10 minutes before using the product again.
- Friction can be harmful to the skin, so do not use it directly. Do not place heavy objects on the product.
- Make sure the product does not come in contact with pointed or sharp objects.

Elektra

- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Do not use the product near a bathtub, shower, swimming pool, or near a tub or laundry tub filled with water. Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get into or on the product, turn the product off immediately and contact our customer service department.
Do not touch the product with wet hands.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Always ensure that the moving parts have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Packaging Contents

- Good4You Infrared Backband - Model MY-BH01

Attached accessories

- 1x User Manual
- 1x USB cable (1m)
- 1x 5000mAh Battery
- 1x Extension Strap


Product Introduction

Infrared heating function


- The waist belt contains an infrared heating system that can cover the entire lower back area and penetrate deep muscles and soft tissues, offering the best heat therapy for your health. If you use the product regularly, your body will experience some pleasant changes, such as improved blood circulation, enhanced muscle flexibility, back pain relief, a cheerful mood, and more. This far-infrared heating belt delivers deeply penetrating heat (infrared rays) that can reach all the way to your bones and create a uniform warming effect without side effects or risks.
This product is also effective for knee relaxation, stomach, and thigh pain relief.
- The product has a built-in 5000mAh lithium battery, making it portable with a long battery life, so you can use it anytime and anywhere.

Product Features


Product Overview

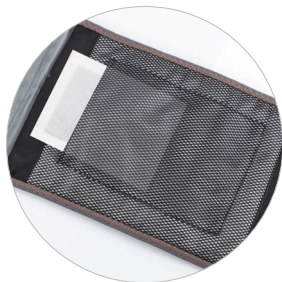Design with mesh pocket

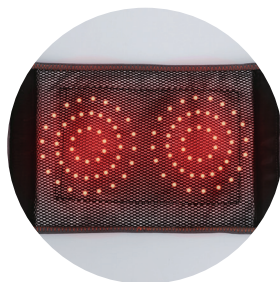Dual LED infrared heating zone

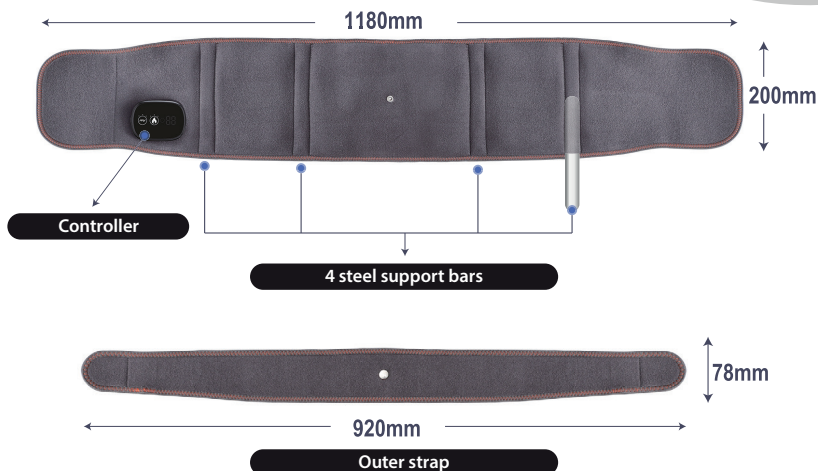Product Features

Battery parameters

Technical Parameters	Type
Charging current	2000mA
Charging time	2-3 hours
Battery capacity	5000mAh


Working parameters

Technical Parameters	Type
Working voltage	5V
Working current	1850mA
Standby current	0.1mA

Battery Life

Technical Parameters	Type
Massage (high -grade)	1440min
Heating (high -grade)	102min
Massage + Heating (high -grade)	66min

Heating module parameters

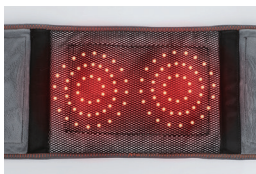
Technical Parameters	Type
Working voltage	5V
Working current	1650mA

Startup instructions

1. Attach the outer belt to the waistband.
 2. Tie the waistband around your waist and adjust the proper tightness according to your waist.
 3. Fasten the outer belt to get a better compression effect.
1. It is recommended to wear the waistband for 30-60 minutes at a time in the beginning. If no side effects or pain are experienced, the usage time can be gradually increased.
 2. Do not wear the waistband while sleeping.
 3. Do not charge it when using it.
 4. It takes a while to heat, please wait patiently.

Instructions for use

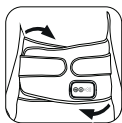01

Tie the waistband around the waist.

02

Long press the power button, the machine beeps and starts up.
If it cannot turn on, try again after charging.

After starting the machine, the default anti-touch mode is disabled and the indicator light is green. You can adjust the massage or heating mode as desired.

03

Press the massage button to turn on the massage function. Press it again to adjust the massage intensity and the corresponding massage indicator will change accordingly.

- The first green indicator light is lit, indicating a low intensity massage.
- The second green indicator light is on, indicating a medium-intensity massage.
- The third green indicator light is on, indicating a high-intensity massage.

No green indicator light means the massage function is disabled.

04

Press the heating button to turn on the heating function. Press it again to adjust the heating mode and the corresponding heating indicator will change accordingly.

- The first green indicator light is on, indicating that the temperature is low.
- The second green indicator light is on, indicating that the temperature is average.
- The third green indicator light is on, indicating that the temperature is high.

No green indicator light means the heating function is turned off.

05

After you set the mode, briefly press the power button. The anti-touch mode turns on and the indicator light turns red. At this time, you cannot freely adjust the massage or heating mode. This prevents accidental touching of the button, which could lead to an inadvertent change of the massage or heating setting during use.

06

In any state, long press and hold the power button to beep and turn the machine off.

Remote Control

Safety Instructions

- When starting to use this product, it is recommended to select the low temperature mode. It is normal for the skin to become slightly red and warm to the touch after use; use the product with confidence. However, if severe redness and swelling occur, discontinue use immediately and consult a physician.
- It is recommended to avoid direct contact with the skin and avoid prolonged use of the high temperature mode to avoid burns.
- It is recommended to start with the lowest temperature and gradually adjust to an appropriate temperature depending on personal circumstances to avoid burning the skin of people who are not sensitive to heat.
- If you experience dizziness, headache, asthma or other discomfort during use, stop using immediately and investigate the cause of the discomfort.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

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