Gebruiksaanwijzing - Roodlichttherapie-apparaat
User Manual - Red Light Therapy Device
Manual d'utilisation - Appareil de luminothérapie rouge
Gebrauchsanweisung - Rotlichttherapiegerät

Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using 'the Red Light Therapy Device', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it!
Read the instructions for use carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- This product is not suitable for individuals (including children) with physical or mental limitations, or for those who lack experience or knowledge of this product. Use of the product is only permitted if supervised by a responsible person who ensures the user's safety. Keep and store this product out of the reach of children. We recommend using this product only if you are 16 years of age or older.
- Do not use this product in combination with implanted electronic devices in the body, such as pacemakers, artificial joints, etc. Consult your doctor if in doubt.
- Do not use this product if you are pregnant, have a severe heart condition, are seriously ill, or are an emergency patient. Consult your doctor if unsure.
- Do not use this product near babies.
- Do not use this product while sleeping.
- Do not use this device for more than 20 minutes in one therapy session.
- Medicines, cosmetics, or dietary supplements may cause hypersensitivity or allergic reactions on the skin. Stop the treatment immediately if this occurs.
- Do not use this product on open or fresh wounds.
- Do not use this product in combination with ointments, balms, or creams that produce heat. This may cause skin burns.
- If you lend this product to someone else, make sure to provide the instruction manual with it.

Elektra

- Do not disassemble, repair, or modify this product yourself. Doing so may damage the product. Warranty will be void if the product is opened or repaired by unauthorized persons. For your safety, such repairs must only be performed by authorized service staff. Always return the product to customer service to avoid damage or personal injury.
- Do not use this product near a bath, shower, swimming pool, or any other container filled with water.
- Ensure this product does not come into contact with water. Do not immerse the product in water. If water does get into the product, stop using it immediately and contact customer service.
- Do not touch the product with wet hands.
- Do not cover the product while it is in use. Never place pillows, blankets, or other objects on top of the device as this may cause fire, electric shock, personal injury, or damage.
- Ensure ventilation openings remain unobstructed during use.

Other

- Ensure that the buttons are not pressed unintentionally when plugging in or moving the product.
- Do not insert objects into the openings or slots of the product. Ensure that no foreign objects enter the product housing.

About the Red Light Therapy Device

The red light therapy device uses a combination of red light (660 nm) and near-infrared light (850 nm). These specific wavelengths are considered ideal for supporting recovery, improving skin health, and promoting overall well-being. The device contains 80 lamps with visible red light (660 nm) and 80 lamps with invisible near-infrared light (850 nm).

Light therapy lamps are designed to target the mitochondria - the energy-producing units of our cells - to stimulate them. By improving mitochondrial function, the body's overall health can be enhanced. The LEDs are tuned to therapeutic wavelengths (expressed in nanometers) to ensure maximum absorption of light in the body. The lamp emits two different types of light: red light and near-infrared light. Research has shown that this combination is highly effective, as each wavelength has its own working depth. Red light penetrates the skin up to a certain depth, while near-infrared light penetrates deeper into muscles, joints, and tissues. Near-infrared light (short wavelength) is hardly visible but penetrates deepest and effectively reaches muscles, joints, and deeper tissues.


Benefits

RED - Red Light 660 nm (visible)

- Penetrates deeply into the skin
- Relaxes muscles
- Firms the skin
- Reduces wrinkles and age spots
- Improves blood circulation and increases oxygen uptake in the blood
- Reduces skin problems and conditions
- Reduces fatigue symptoms
- Supports sports performance
- Promotes better sleep quality
- Relieves stress

NIR - Near-Infrared Light 850 nm (invisible)

- Penetrates deeply into muscles, joints, and tissues
- Reaches IR-A radiation (short wavelength) for maximum penetration
- Relieves pain and reduces inflammation
- Supports weight management
- Stimulates collagen production (anti-aging effect)
- Promotes healing of muscles and tissues
- Reduces symptoms of conditions such as arthritis, fibromyalgia, and ankylosing spondylitis

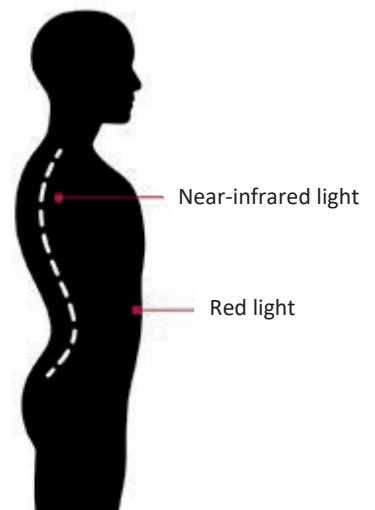Recommended Use

- The ideal usage frequency is between 3 and 7 times per week, maximum once per day per body part.
- Start slowly: position the lamp only for a few minutes at a greater distance so your body can get used to it.
- **Guideline:** the closer the lamp is to the skin, the shorter the recommended treatment duration. Recommended duration is 5 to 20 minutes at a distance of up to 90 cm.
- Example of a recommended treatment: 12 minutes at 30 cm distance.
- Use the lamp for a maximum of 20 minutes. If you want to use the lamp for 20 minutes, increase the distance.
- The lamp has 3 modes: RED infrared light (660 nm), NIR infrared light (850 nm), or a combination of both.
- Note! The lamp **does not** emit heat. The lamp is designed exclusively for LED-light therapy.
- Start the desired program, then put on your protective glasses and try not to look directly into the lamp during treatment.
- Ensure the body part you want to treat is fully exposed to the lamp's light and not covered by clothing.
- **Always consult a doctor before use in case of contraindications, medical conditions, or if in doubt.**


Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents

- Good4You Red Light Therapy Device - Model ZX-303 (A)

Included Accessories

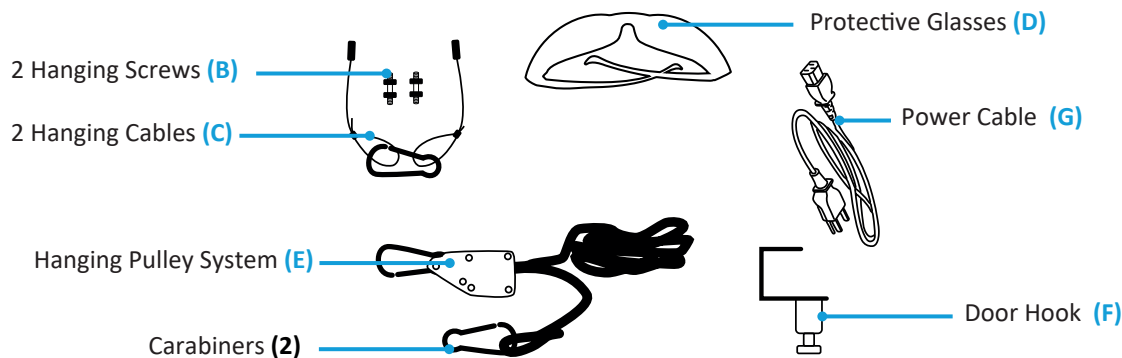Control Panel

Step 1 – Set Session Time

- Press the ▲ button to increase the session time in steps of 5 minutes, up to a maximum of 30 minutes.
- Press the ▼ button to decrease the session time.

Step 2 – Select Program



- After setting the session time, choose the light mode with the “RED” or “NIR” buttons:
 - RED – activates red light (660 nm).
 - NIR – activates near-infrared light (850 nm).
 - Combinatie – activates both red light and near-infrared light simultaneously.
- Refer to the previous page under Benefits for details on the specific advantages of RED light (660 nm) and NIR light (850 nm).

Step 3 – Program Indicators

- Red indicator light (R) = red light is active.
- Green indicator light (IR) = near-infrared light is active.

Note:

Near-infrared light (NIR) is not visible to the naked eye.

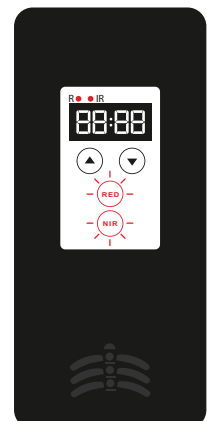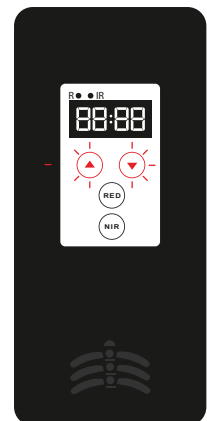Mounting Instructions

i Tip: On the back of the lamp there is a foldable stand.
This allows you to simply place the lamp upright without hanging it.

Do you still want to hang the lamp? Follow the steps below:

Attaching the hanging cables

1. Screw the two hanging screws (B) into the screw holes at the top of the lamp. (Figure 1)
2. Attach the hanging cables (C) to the screws. (Figure 2)

Option A – Hanging on a door hook

1. Place the door hook (F) over the top edge of a standard door.
2. Insert both cable loops through the carabiners and hang these on the door hook. (Figure 3)

Option B – Hanging with pulley system

1. Place the door hook (F) over the top edge of the door.
2. Insert both cable loops through the carabiners of the pulley system (E) and hang the pulley on the door hook. (Figure 4)

Adjusting the height with the pulley system

Raising the lamp

- Support the lamp with one hand and pull the cord with the other until the desired height is reached.
- The pulley will automatically lock. (Figure 5)

Lowering the lamp

- Support the lamp and press the handle to release the cord. (Figure 6)
- Slowly lower the lamp; the pulley will lock automatically.

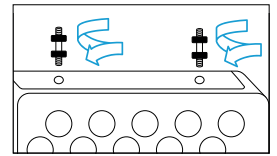
Figure 1

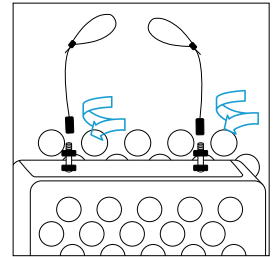
Figure 2

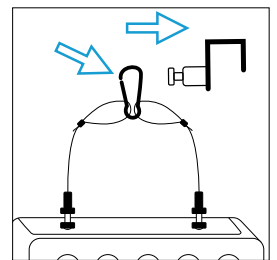
Figure 3

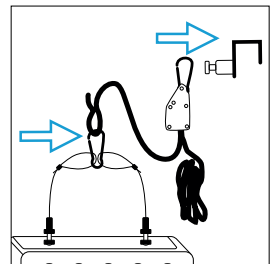
Figure 4

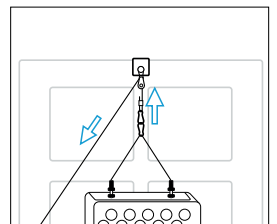
Figure 5

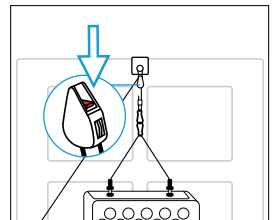
Figure 6

Technical data

Power	400W
Voltage - Frequency	90V-265V - 50/60Hz
Number of LEDs	80 (660nm 80 pcs + 850nm 80 pcs)
Timer settings	5 / 10 / 15 / 20 / 25 / 30 minutes
Product dimensions	36.6 × 21.1 × 6.4 cm
Product weight	2,65 kg

Choosing a treatment option

The following guidelines are general advice for common health and wellness concerns. Your personal treatment needs may vary.

Skin health (wrinkles, acne, wounds)

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 3–5 times
- **Program:** Red light only (RED button)
- **Target area:** All skin areas with visible complaints
- **Distance from the Good4You lamp:** 10 cm

Joint and arthritis pain

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 6–7 times (if severe, possibly twice daily)
- **Program:** Near-infrared light only (NIR button)
- **Target area:** Painful joint
- **Distance from the Good4You lamp:** 10 cm

Weight loss

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 5–7 times (best results before or after exercise)
- **Program:** Red light + near-infrared light (RED & NIR button)
- **Target area:** Area with cellulite or fat, or the entire body
- **Distance from the Good4You lamp:** 10 cm

Muscle recovery

- **Treatment time per session:** 10 minutes
- **Sessions per week:** 3–5 times
- **Program:** Red light + near-infrared light (RED & NIR button)
- **Target area:** Stressed muscles (e.g., after sports)
- **Distance from the Good4You lamp:** 10 cm

Sleep and circadian rhythm

- **Treatment time per session:** 10–15 minutes
- **Sessions per week:** 5–7 times (preferably 30 minutes before sleeping)
- **Program:** Red light only (RED button)
- **Target area:** Whole body; place the lamp where you sleep
- **Distance from the Good4You lamp:** Free to determine in the bedroom

Brain health and cognitive support

- **Treatment time per session:** 10 minutes, twice daily
- **Sessions per week:** 5–7 times
- **Program:** Near-infrared light only (NIR button)
- **Target area:** Head
- **Distance from the Good4You lamp:** 10 cm

Boosting Testosterone

- **Treatment time per session:** 5–10 minutes, 2× per day
- **Sessions per week:** 5–7 times
- **Program:** Near-infrared light only (NIR button)
- **Target area:** Testicles
- **Distance from the Good4You lamp:** 10 cm

Thyroid Function

- **Treatment time per session:** 8–10 minutes
- **Sessions per week:** 5–7 times
- **Program:** Near-infrared light only (NIR button) or red + near-infrared light (RED & NIR button)
- **Target area:** Thyroid area (slightly tilt head backward)
- **Distance from the Good4You lamp:** 10 cm

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.

- Turn off and unplug the product before cleaning. Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Extend warranty for FREE

1) Scan the QR code below to find Good4You's Facebook page and then 'Like' it.

2) Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 1 year to 3 years.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

Troubleshooting

Problem	Possible Causes	Solution
Lamp does not work	1 The lamp is not switched on. 2 The adapter is not properly connected.	1 Press the red ON/OFF switch at the back of the lamp to switch on the LED lights and start the session. 2 Check if the adapter is properly connected to both the lamp and the power socket.
Lamp suddenly stops working	The lamp automatically switches off after the set working time. If no time is set, the lamp switches off after the default working time of 10 minutes.	Switch the lamp back on to start a new red light therapy session.
LED lights are not on while the lamp is working	Near-infrared light (NIR) is not visible to the naked eye.	Check if the program indicator lights above the display are on for the selected program.

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© OnlineGoods