
Gebruiksaanwijzing - Torch Roodlichttherapie
User Manual - Torch Red Light Therapy
Manual d'utilisation - Lampe de thérapie par lumière rouge
Gebrauchsanweisung - Taschenlampe Rotlichttherapie

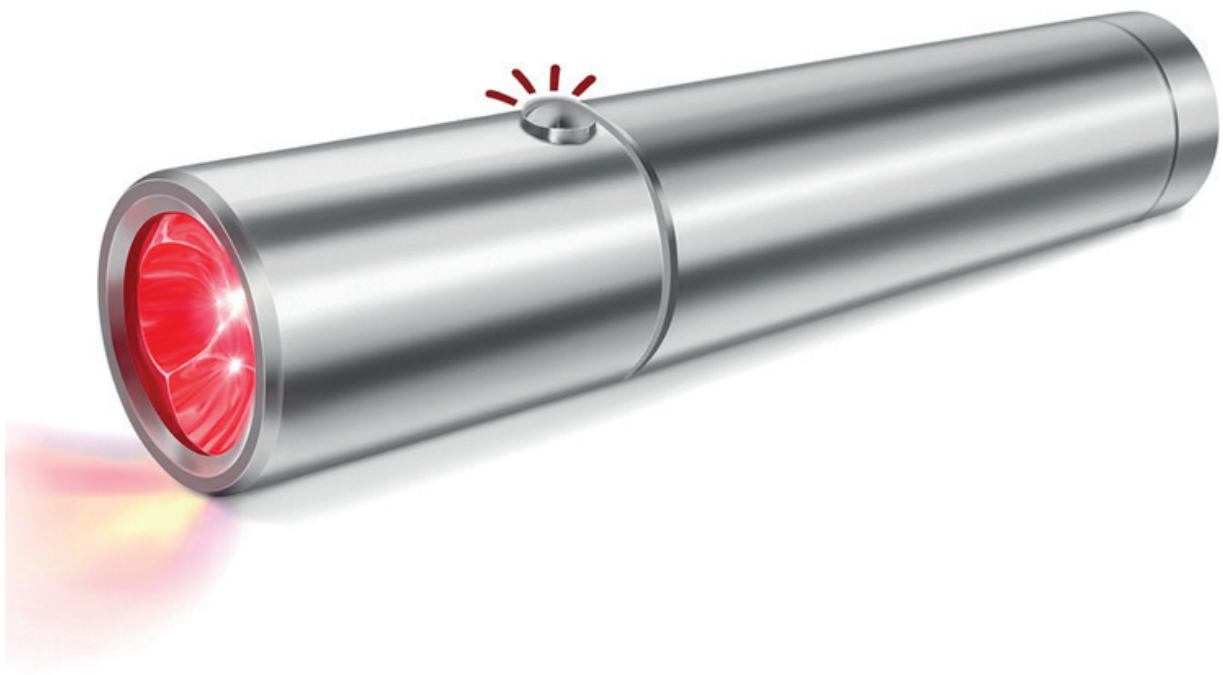Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using 'the Torch Red Light Therapy', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it!
Read the instructions for use carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or those lacking experience or knowledge of the product. Use is only permitted under supervision of a responsible person ensuring user safety. Keep out of reach of children. Recommended for use only by individuals aged 16 or older.
- Do not use the product if you have implanted electronic devices, such as an artificial heart or lungs. Consult your doctor if in doubt.
- Do not use the device if pregnant, with severe heart disease, in serious illness, or in emergency conditions. Consult your doctor if in doubt.
- Do not use the product near babies.
- Do not use the product for more than 20 minutes in one therapy session.
- Medication, cosmetics, or foods may cause hypersensitivity or allergic skin reactions under certain conditions. Stop using the device if this occurs.
- Do not use the product on open or fresh wounds.
- Do not use the product with ointments, creams, or balms containing heat-producing ingredients, as this may cause skin burns.
- If lending the product to someone else, always provide the user manual along with it.

Elektra

- Do not disassemble, repair, or modify the product. Doing so may cause damage. Warranty is void if you open or repair the product yourself. For safety, parts may only be replaced by authorized personnel. Always return the product to a service center for repairs to prevent damage or personal injury.
- Do not use the product near a bath, shower, swimming pool, or sink filled with water.
- Ensure the product does not come into contact with water. Do not immerse in water. If water does enter the product, switch it off immediately and contact Customer Service.
- Do not touch the product with wet hands.

Other

- Do not insert objects into the openings or slots of the product. Ensure that no foreign objects enter the product housing.

About the Torch

The Torch uses a combination of red light (630 nm + 660 nm) and near-infrared light (850 nm). These specific wavelengths are considered ideal for supporting recovery, skin improvement, and overall health.

Light therapy lamps are designed to stimulate the mitochondria – the energy-producing units of our cells. By improving mitochondrial function, the overall health of the body can be enhanced. Red light reaches depths that UV, blue, and green light cannot, and has a healing effect on the skin. Within the infrared spectrum, the rule applies: the shorter the wavelength, the deeper the radiation penetrates. Near-infrared light (short wave) is barely visible but penetrates the deepest, reaching muscles, joints, and deeper tissues more effectively.

Recommended Use

- Aim the Torch at an uncovered part of your forearm for 5 minutes, at a distance of 2.5 to 15 cm.
- Do not keep the Torch on for more than 5 minutes at a time. If you want to treat multiple areas of the body, switch the device off for at least 1 minute after each 5-minute session. This prevents overheating.
- The treatment area is the part illuminated by the red light beam. A small joint, such as a finger, only needs to be treated once, as the light covers the entire joint.
- Light therapy works best when applied consistently. Best results are achieved with daily use, but it is recommended to use the device at least 4 times per week. For severe injuries, chronic pain, or serious skin conditions, using the Torch up to 3 times per day may be beneficial.
- Most users notice benefits within a week of consistent use. For some, it may take up to 3 months before results become visible. Daily use maximizes the benefits of the therapy.
- Note! The Torch **does not** emit heat. The lamp is designed solely for LED light therapy.
- Ensure that the area you want to treat is fully exposed to the Torch light and not covered by clothing.
- **Always consult a doctor before use in case of contraindications, medical conditions, or doubt.**

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents

- Good4You Torch Red Light Therapy - Model AL09

Included Accessories

- 1x User Manual
- 1x Protective Glasses
- 1x Charger
- 2x Rechargeable batteries
- 1x Travel case

Startup Instructions


Step 1:

Connect the charger to your laptop or a suitable power adapter. Then place the battery in the charger and charge it for at least 1 hour before using the device.


Step 2:

Unscrew the housing of the Torch, insert the charged battery with the positive pole (+) facing up and the negative pole (–) facing down, and screw the housing back on tightly.

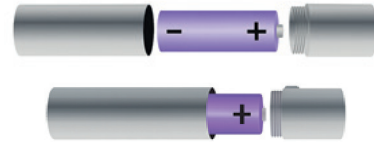
Step 3:

Press the silver button on the Torch and hold it about 2.5 cm (1 inch) from the treatment area (bare skin) for 5 minutes. (The Torch switches off automatically after 5 minutes.)

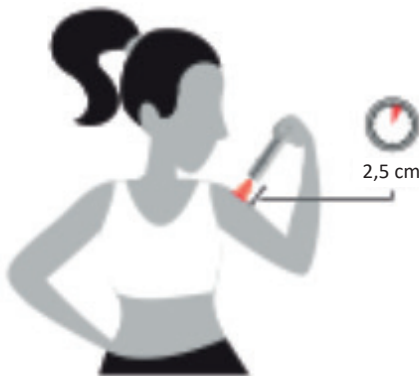Step 4:

Press the silver button again to end the session if the red light is still on. If the session ends automatically, no further action is required.

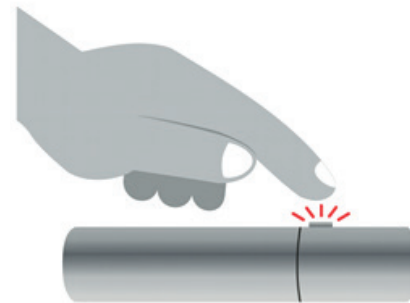Technical data

Power	9W
Irradiance	>133 mw/cm2
Product dimensions	2.5 x 12.5 cm x 0.5 cm
Product weight	200 g

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.

- Turn off and unplug the product before cleaning. Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Extend warranty for FREE

1) Scan the QR code below to find Good4You's Facebook page and then 'Like' it.

2) Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 1 year to 3 years.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

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© OnlineGoods