




Air Compression Leg Massager

User Guide

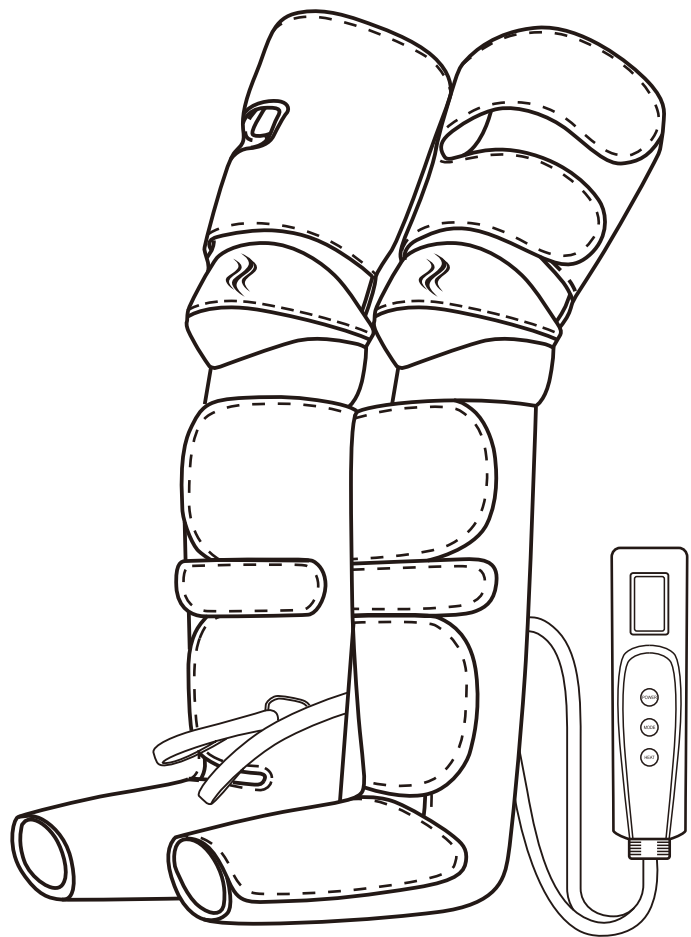

NEED HELP?

WE'RE HERE TO HELP!

Thank you for purchasing the Air Compression Massager. This User Guide is intended to provide you with guidelines to ensure that operation of this product is safe and does not pose risk to the user. Any use that does not conform to the guidelines described in this User Guide may void the limited warranty.


Please read all directions before using the product and retain this guide for reference. This product is intended for household use only. This product is covered by a limited one-year warranty. Coverage is subject to limits and exclusions. See warranty for details.

Box Contents

- Leg massager
- Adapter
- User guide

Parts

1. Leg massager
2. Air pressure hose
3. LCD Screen
4. Control panel

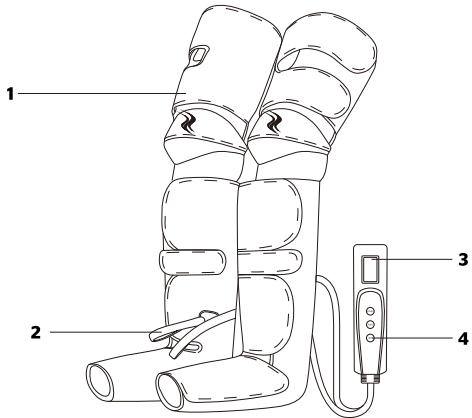


Safety Precautions

- Do not use voltage beyond the requirement of this product.
- Do not attempt to repair or dismantle the product. There are no user-serviceable parts.
- Do not use the product in a place where it may become wet or damp. If the unit becomes wet, do not use.
- Stop using the product immediately if you feel uncomfortable or any abnormalities occur.
- Do not use this product in a high-temperature environment.
- Do not use this product if you are currently pregnant or menstruating.
- Do not place the unit in a place where it may be in contact with flammable items, explosive items, or chemicals.
- Do not touch the unit with any sharp object, or risk damage.
- Do not use this product if any of the below apply to you:
 - Skin is infected, swollen, burnt, or otherwise damaged.
 - You have a pacemaker or other medical implant.
 - You are suffering from heart disease.
 - Your bones are fragile or injured.
 - You are suffering from cancer.
 - You are receiving treatment or your blood pressure is not normal.
 - You are suffering from osteoporosis.

Controls

1. <POWER> - Press 3s to open the machine, Press second time to adjust intensity to LOW-MIDDLE-HIGH, Press 3s again to off the product
2. <MODE> - Have 3 mode can be choose. Mode 1: Calf + Foot Airbag, Model 2: Only Calf, Mode 3: Only Foot.
3. <HEAT> - Heat ON/OFF, the heat place is in the bottom of knee.

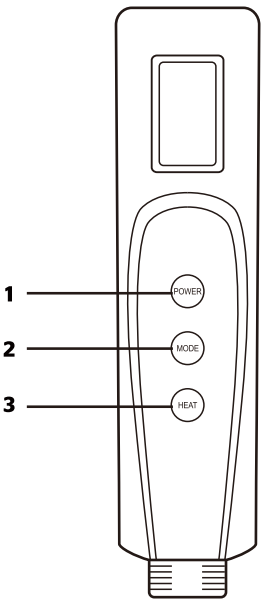

Operation

1. Slip on each massager, keeping your feet flat on the bottom. Note that each is marked L and R for each foot.
2. Use the strap to tightly secure your foot into the massager.
3. Use the strap to tightly secure your leg into the massager.
4. Move your foot to make sure that the massager is securely fastened.
5. Insert the two air pressure connections into the controller.
6. Press the power button While on the LED will stay on
7. Adjust the air pressure strength and mode by your hand, and can see your operation in LCD screen

NOTE: the device will automatically shut off after 20 minutes.

For best results, do not use the leg massager for more than an hour.

Care and maintenance

- Store the product in a dry place, away from moisture and direct sunlight.
- Wipe clean with a dry towel.

Troubleshooting

Problem	Cause	Solution
The device does not work.	The unit is not connected to a power source.	Make sure the unit is fully connected to the power.
	The power button was not pressed.	Press the power button.
The device suddenly stopped working.	The unit was running more than 15 minutes.	Restart the unit to resume use.

Specifications

Rated voltage: 12V 2A
Power consumption: 24W
Default time: 20 min.