
Gebruiksaanwijzing - Shiatsu Voetmassageapparaat
User Manual - Shiatsu Foot Massager
Manual d'utilisation - Shiatsu Masseur de pieds
Gebrauchsanweisung - Shiatsu-Fußmassagegerät

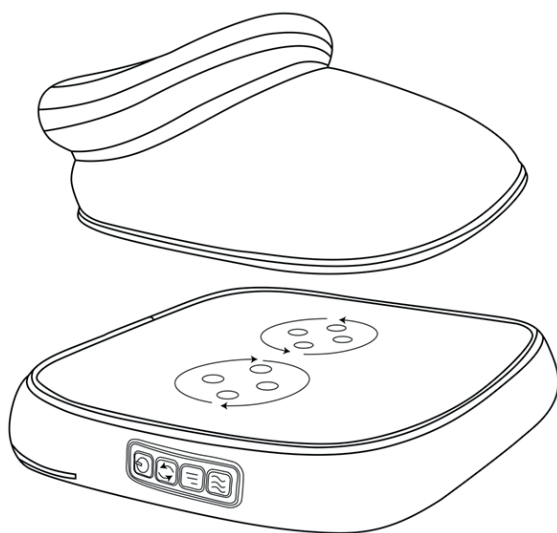Language:    

Model: FY-5007

Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl

Intended Use

- This product is intended only for the purposes described in the instructions for use.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using "the Foot Massager", hereinafter referred to as "the product".

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They may choke on it! Read the instructions for use carefully before use and keep them carefully!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Be aware of the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or people sen who have no experience with or knowledge about the product. Use of the product is permitted only when a supervisor is present and responsible for the safety of the user. Use and store the product out of the reach of children. We recommend using the product only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; are pregnant; have injuries or lesions to your foot(s), or have conditions such as warts, swimmer's eczema, rashes or open wounds; have a malignant tumor, heart disease, osteoporosis or a broken spine; Have had surgery on your foot(s); Have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nervous disorders, acute inflammation, skin disease, broken skin, or pain with no obvious or immediate cause.
- Do not use the product with more than one person at a time. Do not use the product near your eyes or other sensitive areas. Do not use the product when you are in bed or sleeping. Do not use the product just before you go to bed; the product has a stimulating effect. Do not use the product when you are driving a vehicle or operating machinery. Do not use the product vertically to avoid accidents or injuries to the turbine, malfunctions, etc.
- If you lend the product to someone else, be sure to give the product with the manual together.
- Use the product on a level surface. If it is not stable, the body will not be balanced and an accident can be caused.
- The product is special for the legs, feet, etc., do not use the waist, neck and other parts outside the purpose of use.

- Prolonged use of the product may cause overheating. After 15 minutes of use, insert a break of at least 15 minutes before using the product again.
- Do not leave the product running unloaded or reverse the direction of rotation of the product.
- Friction can be harmful to the skin, so do not use it directly.
- Do not place heavy objects on the product.
- Make sure the product does not come in contact with pointed or sharp objects.

Medical Notes

- The product is not a substitute for treatments from a physician or physical therapist.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment or are using medical equipment.
- To avoid excessive over stimulation of the muscles and nerves, it is recommended that you do not use the product for more than 15 minutes at a time.
- If you experience pain or discomfort while using the product or find the massage unpleasant, discontinue use and consult your physician.

Electra

- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Do not use the product near a bathtub, shower, swimming pool, or near a tub or laundry tub filled with water.
- Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get into or on the product, turn the product off immediately and contact our customer service department.
- Do not touch the product with wet hands.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Always ensure that the moving parts have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents

- Good4You Shiatsu Foot massager - Model FY-5007

Attached accessories


- 1x User Manual
- 1x Power Adapter

Instructions for use

- ⚠ Prolonged use of the product may lead to overheating. After 15 minutes of use, take a break of at least 15 minutes before using the product again and do not use it more than 3 times a day.
- ⚠ Do not cover the product when it is on. Never place pillows or blankets on or against the unit. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- ⚠ Keep the product away from moisture, heat, areas where the temperature is relatively high, and heat sources such as ovens.

Start-up instructions

1. Connect the Dc-Connector to the device plug.
2. Plug the connector into the power outlet.
3. Press the "⏻" button to turn on the massage cushion.
4. Press the "↺" button to adjust the massage direction.
5. Press the "|||" button to adjust the intensity. The massage speed goes from low to medium and then to high.
6. Press the "⌘" button to turn the heat mode on/off.
7. Press the "⏻" button to turn off the massage cushion. Unplug it after use.

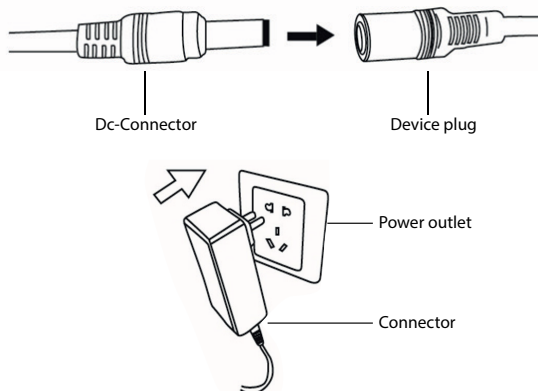Operation panel

On/Off: Starts or ends the massage. After 15 minutes, it stops automatically.

Massage direction: Choose whether the massage heads rotate counterclockwise or

Intensity: Adjust the rotation speed to your personal preference.

Heat mode: Turns the heat mode on and off.

Technical data

Dimensions	32 x 33 x 19 cm
Weight	1.5 kg
Power	24W
Power - Frequency - Output	12V - 50Hz - 2A
Timer	15 min

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

- Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.
- Turn off and unplug the product before cleaning.
- Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.
- The top of the unit has a removable cover that can be washed.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl. You will receive a response within 24 hours.
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be made valid within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Extend warranty for FREE

- 1) Scan the QR code below to find Good4You's Facebook page and then 'Like' it.

- 2) Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 1 year to 3 years.

Troubleshooting

Problem	Solution
A clattering sound is heard in the movement.	This is perfectly normal and is caused by the operation of the product. It is a normal status.
The massage device switches off suddenly	The device automatically turns off after 15 minutes. Continue use by pressing the power button. If the device has been used for 30 minutes, turn the device back on after resting for 15 minutes.
	Overheating protection is activated. The appliance automatically turns off for safety. Let the appliance rest for 15- 30 minutes before using it again.
Het apparaat start niet	Check if the device has been used for more than 15 minutes. If so, turn the unit back on after resting for 15 minutes. In warmer environments, it may take longer for the device to cool down.
	Check if the adapter output is connected to device power inlet.

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© OnlineGoods