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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by the use of the product beyond its intended use.
- Always observe the following precautions when using "the Foot Massager", hereinafter referred to as "the product".

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it! Please read the operating instructions carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Note the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. If you are not using the product, turn the product off.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or person who have no experience with or knowledge about the product. Use of the product is permitted only when a supervisor is present and responsible for the user's safety. Use and store the product out of the reach of children. We recommend using the product only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; are pregnant; have injuries or lesions to your foot(s), or have conditions such as warts, swimmer's eczema, rashes or open wounds; have a malignant tumor, heart disease, osteoporosis or a broken spine; Have had surgery on your foot(s); Have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nerve disorders, acute inflammation, skin disease, broken skin, or pain without obvious or immediate cause.
- Do not use the product with more than one person at a time. Do not use the product near your eyes or other sensitive areas. Do not use the product when you are in bed or sleeping. Do not use the product just before you go to bed; the product has a stimulating effect. Do not use the product when you are driving a vehicle or operating machinery. Do not use the product vertically to avoid accidents or injuries to the turbine, malfunctions, etc.
- If you lend the product to someone else, be sure to give the product with the manual together.
- Use the product on a level surface. If it is not stable, the body will not be balanced and an accident can be caused.
- The product is special for the legs, feet, etc., do not use the waist, neck and other parts outside the purpose of use.
- Prolonged use of the product may lead to overheating. After 10 minutes of use, insert a break of at least 10 minutes before using the product again.
- Do not leave the product running unloaded or reverse the direction of rotation of the product.
- Friction can be harmful to the skin, so do not use it directly.
- Do not place heavy objects on the product.
- Do not allow the product to come into contact with pointed or sharp objects.

Medical Notes

- The product is not a substitute for treatments from a doctor or physical therapist.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment, or are using medical equipment.
- To avoid excessive over stimulation of the muscles and nerves, it is recommended that you do not use the product for more than 10 minutes at a time.
- If you experience pain or discomfort while using the product or find the massage unpleasant, discontinue use and consult your physician.

Elektra

- Do not disassemble, repair or restore the product. This may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Do not use the product near a bathtub, shower, swimming pool, or near a tub or laundry tub filled with water.
- Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- Do not touch the product with wet hands.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Always ensure that the moving parts have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents


- Good4You Foot massager - Model KTR-401

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
- 1x User Manual
- 1x Remote Control
- 3x AA Batteries
- 2x AAA batteries

Product Features

Product overview

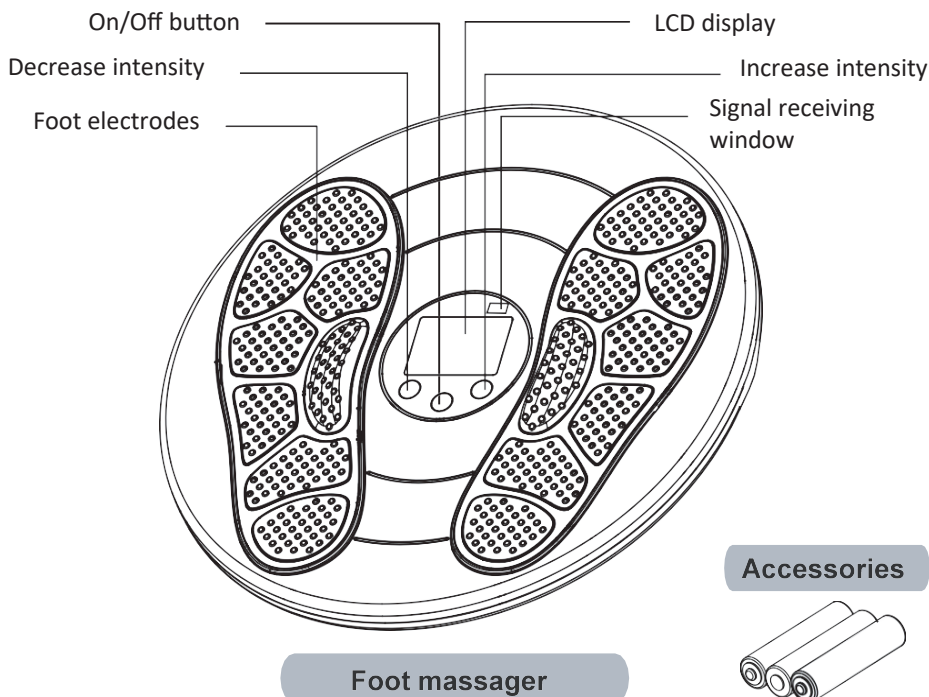
Accessories

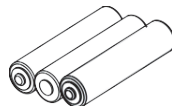
3x AA batteries

Remote control

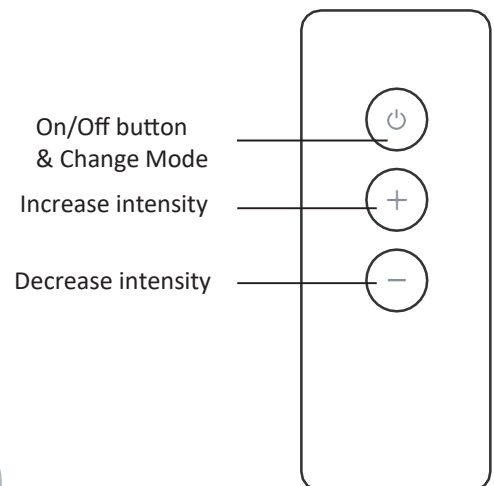
Signal sending window

Accessories

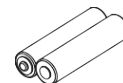2x AAA batteries

LCD display

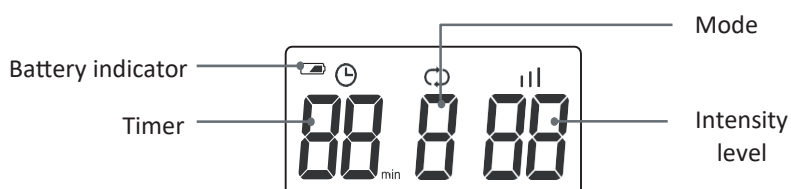Instructions for use

- ⚠ Prolonged use of the foot massager may lead to overheating. After 20 minutes of use, take a break for at least 20 minutes before using it again.
- ⚠ Keep the foot massager and remote control away from moisture, heat, rooms where the temperature is relatively high, and heat sources such as furnaces

4.1 Preparations before use

- 1) Insert 3 AA batteries into the foot massager and 2 AAA batteries into the remote control. Pay attention to the positive and negative poles of the batteries. Close the battery cover after insertion.
- 2) Clean your feet and the foot electrodes with a damp towel and then place your feet on the electrodes; then position your feet so that they make full contact with the two electrodes, then completely relax your body and start operating the remote control. It is recommended to put a flannel on the electrodes in winter to keep your feet warm while using the foot massager.

4.2 Operating the remote control

- 1) Press the power button of the foot massager and the LCD screen will light up, then press the " ⏻ " button of the remote control once to turn on the foot massager, the foot massager will make a "beep" sound and the following default settings will be enabled:

- Timer of 30 minutes
- Mode 1
- Intensity 0

Then you can use the remote control to set the desired mode and intensity; when operating the remote control, point the signal sending window to the signal receiving window on the foot massager.

The foot massager will beep once after successfully receiving the infrared signal. Long press and hold the " ⏻ " button on the remote control to turn off the foot massager.

- 2) To adjust the intensity, press the " + " button to increase the intensity and the " - " button to decrease the intensity. There are a total of 20 different intensity levels.

After turning on the foot massager, it defaults to level " 0 ".

- 3) Press the " ⏻ " button once to switch between 9 different massage modes.

Mode 1 offers an automatic massage, where all 9 different massage modes are combined into 1 massage (everything is massaged), while modes 2 to 9 offer circular massage, i.e. every minute it automatically switches to the next mode; modes 2 to 9 offer single massage where a different point is massaged each time.

*** It is recommended to use the foot massager up to twice a day, 15-30 minutes at a time.**

4.3 Use of batteries

- 1) 3 pieces of AALR6 batteries are required for the foot massager for approximately 60 days of use (based on 30 minutes per day).
- 2) 2 pieces AAA batteries are needed for the remote control for about half a year of use (based on 30 minutes per day).
- 3) If " 🔋 " appears on the LCD display of the foot massager, it means the battery is low. Replace the batteries with new ones.
- 4) If you do not use the foot massager for a long time, remove the batteries, otherwise the device may malfunction due to leaking batteries.
- 5) Do not use batteries other than those specified externally.
- 6) Pay attention to the positive (" + ") and negative (" - ") poles when installing batteries.
- 7) Never use batteries at temperatures above 45°C, otherwise their performance and life will be affected.
- 8) Discarded batteries should be disposed of in accordance with local environmental regulations.

Technical data

Dimensions	36 x 36 x 4 cm
Weight	1,5 kg
Input Massager	DC 12V ~ 1,5A
Input adapter	AC 100-240V ~ 50/60Hz ~ 0,8 A
Number of massage modes	9
Intensity levels	20
Timer	30 min

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

- ⚠ Do not immerse the product in water or other liquids and keep it away from o solvents and harsh cleaning agents.
- Turn off the product before cleaning.
- Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of the product's life, do not throw it away with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects of which the customer was aware at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

Troubleshooting guide

Problem	Cause	Solution
The device does not work	Does the remote control have sufficient battery power?	Replace the batteries with new AAA batteries and insert them with the correct positive and negative poles.
	Does the foot massager have sufficient battery power?	Replace the batteries with new AA batteries and insert them with the correct positive and negative poles.
There is a tingling feeling	Do the electrodes come in contact with your soles or body skin?	The electrodes should make good contact with your soles of your feet or body skin.
	Does the massage take too long?	It is appropriate to massage for 15 to 30 minutes each time, up to twice a day.
	Is your skin too dry?	Wet the soles of your feet with a damp towel.
	Is the intensity of the massage too low?	Increase the intensity of the massage.
	Is the intensity of the massage too high?	Reduce the intensity of the massage.
Low stimulation	Is there dirt or oil stains on the electrodes?	Clean the surface of the hydrogel on the electrodes.
	Is the battery of the foot massager low?	Replace dead batteries with new ones.
Red spots on the skin after use	Is there a strong tingling sensation while using of the product? Did you use an ointment?	Stop using the foot massager immediately and wait for the red spots to disappear after 3-5 days.

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© OnlineGoods