
Gebruiksaanwijzing - Massagemat
User Manual - Massage mat
Manual d'utilisation - Tapis de massage
Gebrauchsanweisung - Massagematte

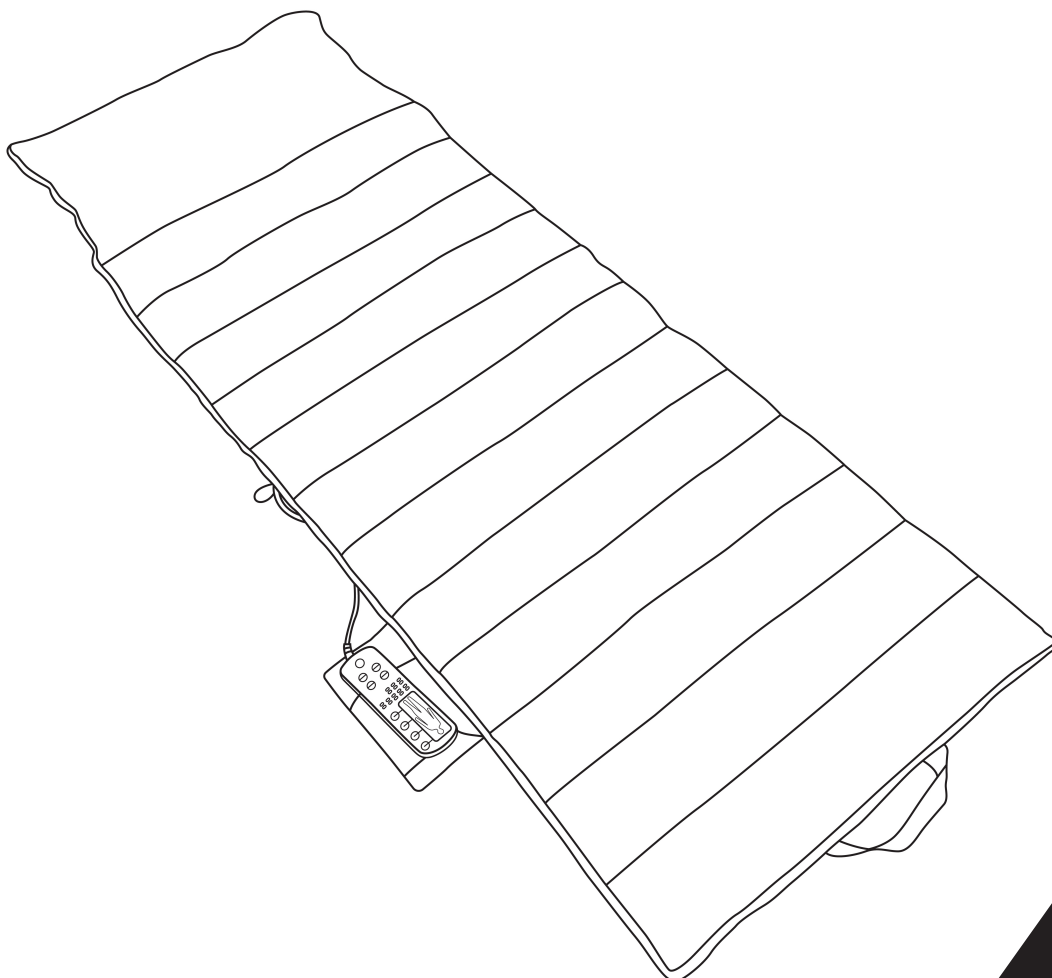Introductie

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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us at email: info@good-4you.nl.

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by the use of the product beyond its intended use.
- Always observe the following precautions when using 'the Massage Mat', hereinafter referred to as 'the Product'.

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it! Read the instructions carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Be aware of the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When you are not using the product, turn the product off and unplug the power cord.
- Do not pull the cord to unplug the product.
- Do not move the product by pulling the cord or twist the cord.
- Make sure you do not trip over the cord and that the cord is not kinked, twisted, or stuck.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or one sen who have no experience with or knowledge about the product. Use of the product is permitted only when a supervisor is present and responsible for the user's safety. Use and store the product out of the reach of children. We recommend using the foot massager only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; are pregnant; have injuries or injuries to your foot(s), or have conditions such as warts, swimmer's eczema, rashes or open wounds; have a malignant tumor, heart disease, osteoporosis or a broken spine; Have had surgery on your foot(s); Have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nervous disorders, acute inflammation, skin disease, broken skin, or pain with no obvious or immediate cause.
- Do not use the product with more than one person at a time.
- Do not use the product near your eyes or other sensitive areas.
- Do not use the product when you are in bed or sleeping.
- Do not use the product just before you go to bed; the product has a stimulating effect.
- Do not use the product when you are driving a vehicle or operating machinery.
- Prolonged use of the product may lead to overheating. After 30 minutes of use, take a break for at least 30 minutes before using the product again.
- Make sure the product does not come into contact with pointed or sharp objects.

Medical Notes

- The product is not a substitute for treatments from a doctor or physical therapist.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment or are using medical equipment.
- To prevent excessive over stimulation of the muscles and nerves, it is recommended that you do not use the product for more than 30 minutes at a time.
- If you experience pain or discomfort while using the product or find the massage unpleasant, discontinue use and consult your physician.

Electra

- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Before connecting the product, ensure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Use only the adapter provided by the supplier.
- Do not use the product near a bathtub, shower, swimming pool, or near a (laundry) tub filled with water.
- Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- Do not touch the product and power cord with wet hands.
- If the power fails during use, turn off the product immediately and unplug the product.
- The product gets hot during use. Do not use the product if you are sensitive to heat. Keep the product away from moisture, heat, areas where the temperature is relatively high, and heat sources such as ovens.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Make sure that the moving parts always have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.


Package Contents

- Good4You Massage Mat - Model SM-603

Attached accessories

- 1x Remote Control
- 1x Power Adapter
- 1x User Manual

Remote control

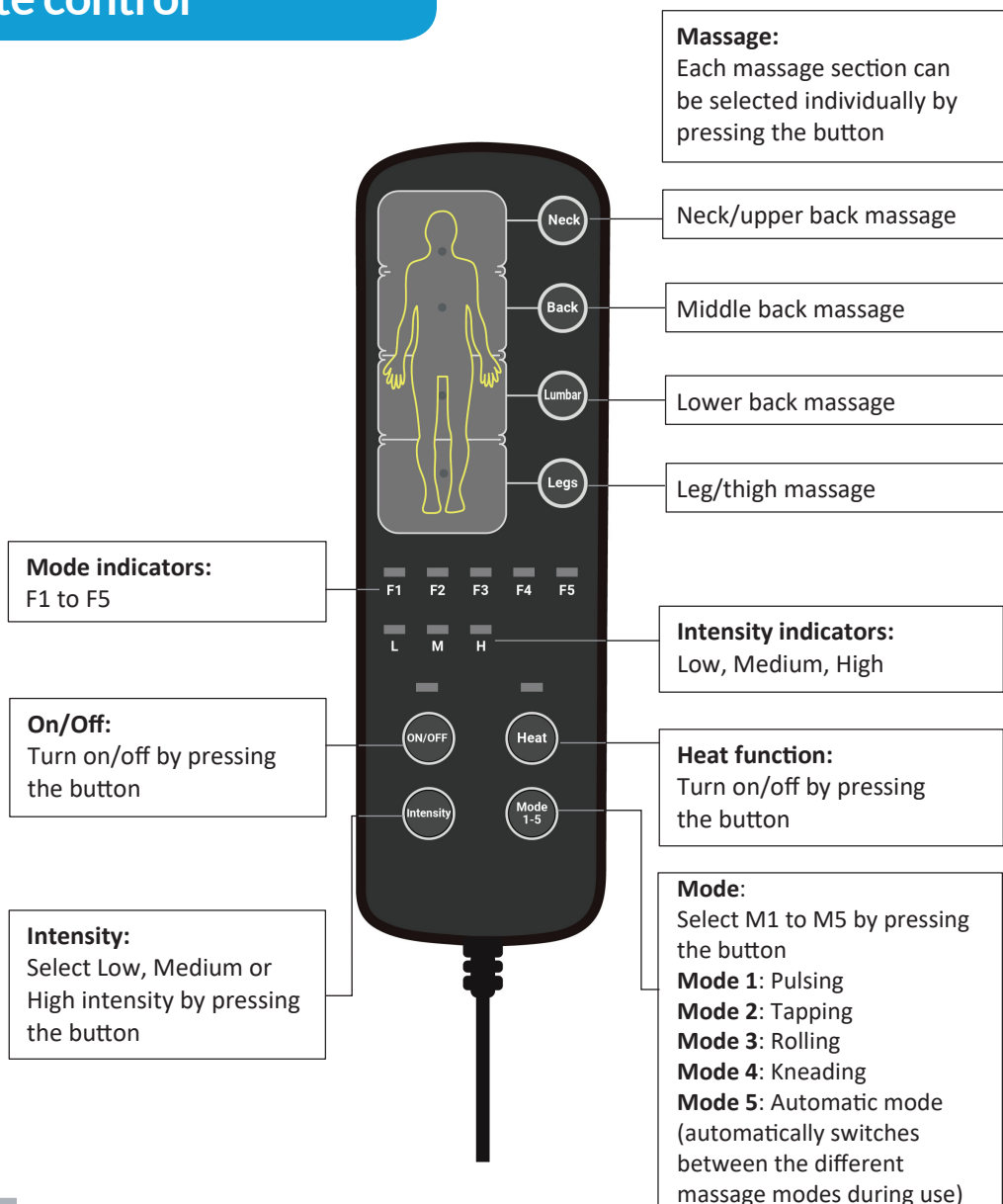
Instructions for use

- ⚠ Prolonged use of the product may lead to overheating. After 30 minutes of use, take a break of at least 30 minutes before using the product again.
- ⚠ Do not cover the product when it is on. Never place pillows or blankets on or against the unit. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- ⚠ Keep the product away from moisture, heat, areas where the temperature is relatively high, and heat sources such as furnaces.

Startup instructions

1. Place the massage mat on a flat and stable surface.
2. Connect the Dc-Connector to the device plug.
3. Plug the plug into the power outlet.
4. The F1 and M LED indicators turn on by default after power on and will turn orange.
5. Select from the 4 massage sections by pressing the corresponding button. Adjust the section by selecting each combination or all 4 at once. The corresponding LED indicator will turn green.
6. Press the "Mode" button to change the mode. The corresponding LED indicator will turn orange.
7. Press the "Intensity" button to change the intensity. The corresponding LED indicator will turn orange.
8. Press the "Heat" button to turn on the heat function. The corresponding LED indicator will turn red. The massage mat is protected against overheating and will automatically turn off if the temperature is too high.
9. Press the "On/Off" button to turn off the massage mat. Unplug it after use and place the remote control in the pocket on the side of the massage mat.

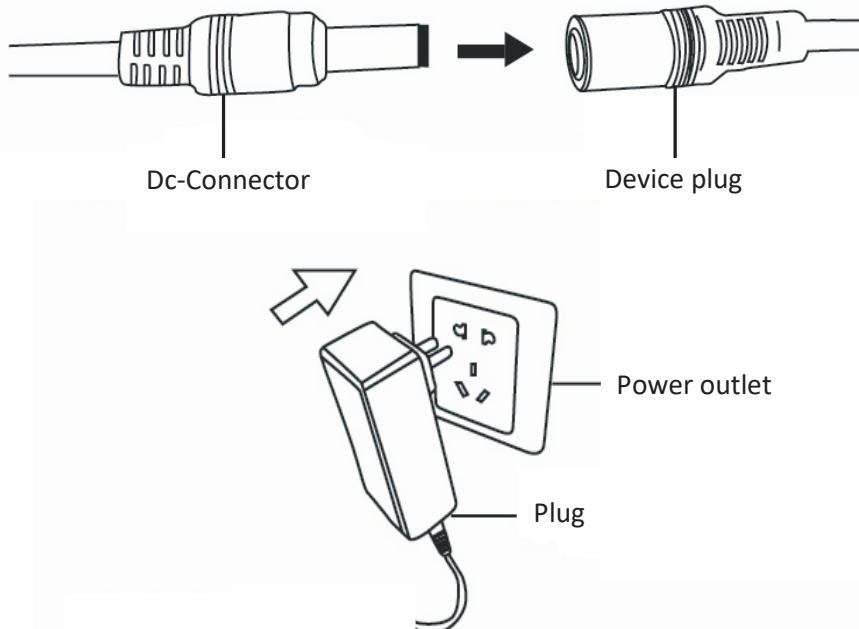Technical data

Dimensions	170 x 55 x 4 cm
Weight	2,8 kg
Power Output	30W
Power supply - Frequency - Output	12V - 50Hz - 2.5A
Number of massage motors	10
Number of intensity indicators	3
Number of massage sections	4
Number of mode indicators	5
Timer	30 min

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

- ⚠ Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents. Turn off and unplug the product before cleaning.
- Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must be able to show proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

Troubleshooting

Problem	Solution
The massage is not active in all massage sections	Possibly the massage motor for the massage section is turned off. Press the button associated with the section to activate it.
The device turns off suddenly	The device automatically turns off when the timer has expired. Continue use by pressing the power button. If the unit has been used for 60 minutes, turn the unit back on after resting for 15 minutes.
	Overheating protection is activated. The device automatically shuts off for safety. Let the device rest for 15-30 minutes before using it again.
The device does not start	Check if the device has been used for more than 30 minutes. If so, turn the device on after resting for 15 minutes. In warmer environments, it may take longer for the device to cool down.
	Make sure the adapter is plugged in.
	Make sure the device plug is connected to the DC connector.
Jerky movement due to vibration of the device	Too much pressure may have been applied to the device. Reduce the load and try again.

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© OnlineGoods