
Gebruiksaanwijzing - Voetmassageapparaat
User Manual - Foot massager
Manual d'utilisation - Appareil de massage des pieds
Gebrauchsanweisung - Fußmassagegerät

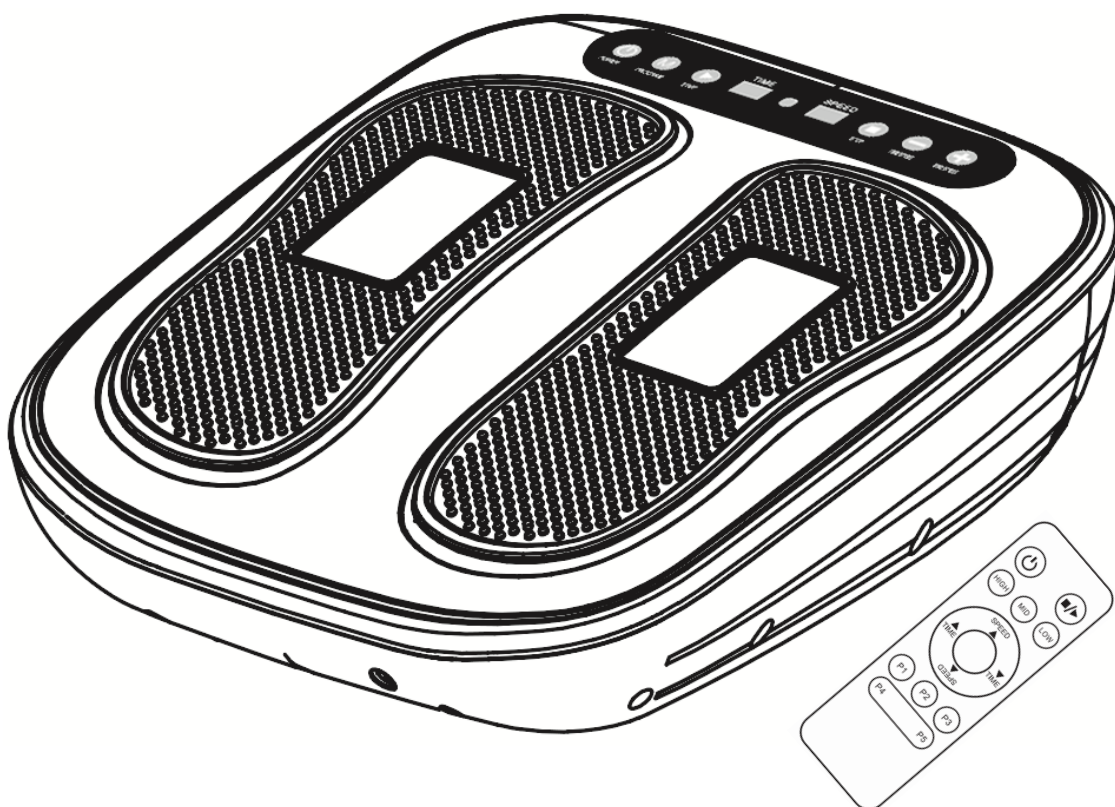Introductie

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Introduction

Thank you for choosing our product! Please always follow the safety instructions for safe use. Should there be any questions about the product please feel free to contact us via email: info@good-4you.nl

Intended Use

- This product is intended only for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by the use of the product beyond its intended use.
- Always observe the following precautions when using "the Foot Massager", hereinafter referred to as "the product".

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They could choke on it! Read the instructions carefully before use and keep them in a safe place!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Be aware of the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When you are not using the product, turn the product off and unplug the power cord.
- Do not pull the cord to unplug the product.
- Do not move the product by pulling the cord or twist the cord.
- Make sure you do not trip over the cord and that the cord is not kinked, twisted, or stuck.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or People who have no experience with or knowledge about the product. Use of the product is permitted only when a supervisor is present and responsible for the user's safety. Use and store the product out of the reach of children. We recommend using the foot massager only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; are pregnant; have injuries or injuries to your foot(s), or have conditions such as warts, swimmer's eczema, rashes or open wounds; have a malignant tumor, heart disease, osteoporosis or a broken spine; Have had surgery on your foot(s); Have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nervous disorders, acute inflammation, skin disease, broken skin, or pain with no obvious or immediate cause.
- Do not use the product with more than one person at a time. Do not use the product near your eyes or other sensitive areas. Do not use the product when you are in bed or sleeping. Do not use the product just before you go to bed; the product has a stimulating effect. Do not use the product when you are driving a vehicle or operating machinery. Do not use the product vertically to avoid accidents or injuries to the turbine, malfunctions, etc.
- If you lend the product to someone else, be sure to give the product with the manual together.
- Use the product on a level surface. If it is not stable, the body will not be balanced and an accident can be caused.

- The product is special for the legs, feet, etc., do not use the waist, neck and other parts outside the purpose of use.
- Prolonged use of the product may cause overheating. After 10 minutes of use, insert a break of at least 10 minutes before using the product again.
- Do not leave the product running unloaded or reverse the direction of rotation of the product.
- Friction can be harmful to the skin, so do not use it directly.
- Do not place heavy objects on the product.
- Make sure the product does not come in contact with pointed or sharp objects.

Medical Notes

- The product is not a substitute for treatments from a doctor or physical therapist.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment or are using medical equipment.
- To prevent excessive over stimulation of the muscles and nerves, it is recommended that you do not use the product for more than 10 minutes at a time.
- If you experience pain or discomfort while using the product or find the massage unpleasant, discontinue use and consult your physician.

Electra

- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Before connecting the product, ensure that the mains voltage indicated on the rating plate matches that of your mains supply.
- Check both power cord and plug regularly for signs of damage or wear, and if there is any, do not use the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- Use only the adapter provided by the supplier.
- Do not use the product near a bathtub, shower, swimming pool, or near a (laundry) tub filled with water.
- Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- Do not touch the product and power cord with wet hands.
- If the power fails during use, immediately turn off the product and unplug the power cord.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Always ensure that the moving parts have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.


Contents of Package


- Good4You Foot massager - Model ETF005

Attached Accessories

- 1x User Manual
- 1x Power Adapter
- 1x Remote Control

Product Features

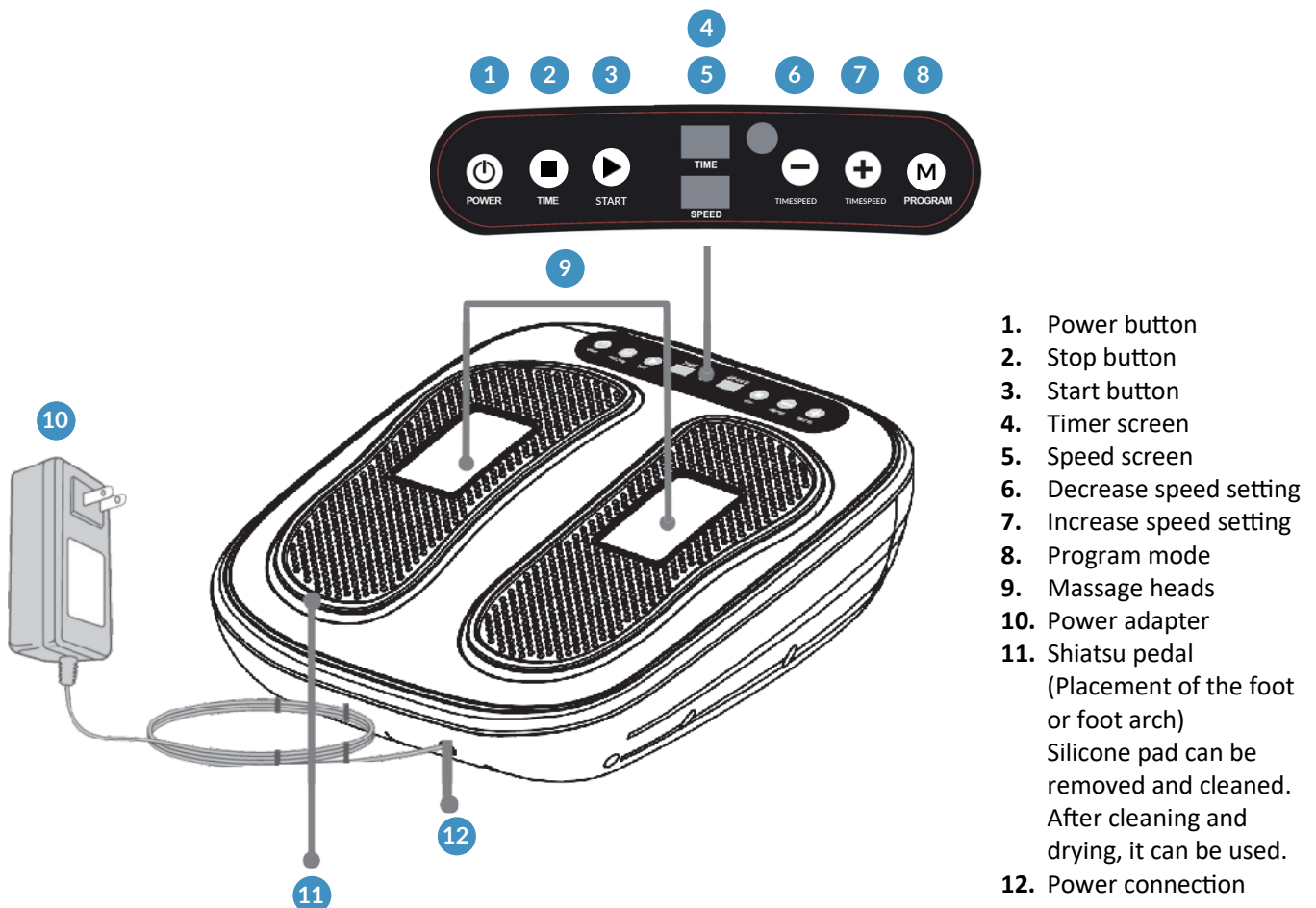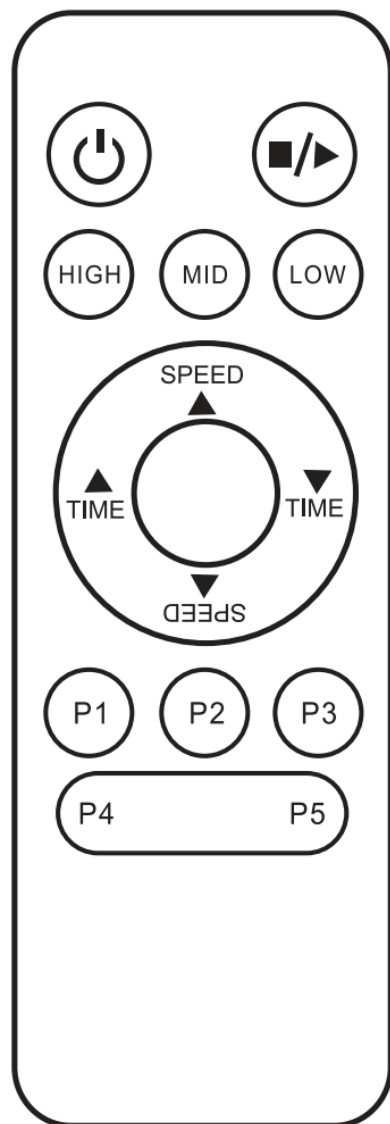 On/Off button

 Start/stop-button

HIGH Select 16 speeds

MID Select 10 speeds

LOW Select 6 speeds

→ These three functions for speed are only available in P1 mode

SPEED ▲ Increase speed

SPEED ▼ Decrease speed

→ The speed range is adjustable from 1 to 20
 Adjustable in P1 mode only

TIMER ▲ Increase time of timer

TIMER ▼ Decrease time of timer

→ The time range is adjustable from 1 to 10 minutes and can only be adjusted after pressing the Start/Stop button
 Adjustable in P1 mode only
 After power on, the default time is 10 minutes

P1

Default time: 10 minutes and adjustable by pressing the Start/Stop button

Function: Vibration only

Speed: Adjustable by pressing Start/stop button

* P1 mode is turned on by default after power on

P2

Default time: 10 minutes

Function: Vibration only

Speed: Speed changes every second for 60 seconds total, from 1 to 20 and 20 to 1

P3

Default time: 10 minutes

Function: Vibration with massage

Speed: The speed changes every second for 60 seconds in total, from 1 to 20 and 20 to 1

P4

Standard time: 10 minutes

Function: Vibration with or without massage

Speed: The speed changes every second for 60 seconds in total, from 1 to 20 and 20 to 1

1st minute: Vibrations only

2nd minute: Vibrations with massage

3rd minute: Only vibrations

4th minute: Vibrations with massage

5th minute: Only vibrations

6th minute: Vibrations with massage

7th minute: Only vibrations

8th minute: Vibrations with massage

9th minute: Only vibrations

10th minute: Vibrations with massage

P5

Standard time: 10 minutes

Function: Vibration with or without massage

Speed: Speed changes every second for 60 seconds total, from 1 to 20 and 20 to 1

2 minutes: Vibrations only


3 minutes: Vibrations with massage

2 minutes: Vibrations only


3 minutes: Vibrations with massage

Suitable for

Feet

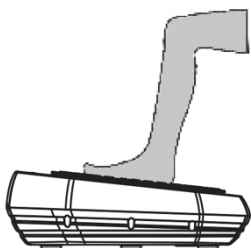Calves

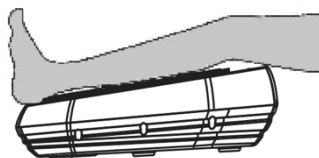Back

⚠ **Note:** always use the product while sitting, never standing!

Technical data

Dimensions	40 x 38 x 12 cm
Weight	4,1 kg
Color	Black
Material	ABS
Cable length	1,8 m
Input Massager	DC 12V2A
Power	30W
Input power adapter	AC 100-240 V ~ 50/60Hz - 0,7 A
Number of positions	5
Speed levels	20
Automatic running time	10 min
Massage type	Shiatsu massage
Suitable for (body part)	Feet, calves and back
Suitable for (functions)	Muscle pain, stiff muscles, stimulate blood circulation, make muscles supple and flexible

Instructions for use

- ⚠ Prolonged use of the product may lead to overheating. After 10 minutes of use, take a break for at least 10 minutes before using the product again.
- ⚠ Do not cover the product when it is on. Never place pillows or blankets on or against the unit. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- ⚠ Keep the product away from moisture, heat, areas where the temperature is relatively high, and heat sources such as furnaces.

1. Device:

1.1 When the plug is plugged in, then press the "POWER"-button and the display will show:

"Default time: 10 minutes. Speed: 00."

1.2 Press the "PROGRAM"-button to view the modes P1, P2, P3, P4 and P5.

Select the desired mode and click "START" (program P2,- P3,- P4- and P5 are selectable with the remote control).

1.3 If you do not want to follow the "1.2" mode, press the "START"-button, the device will start from P1 mode by default. The speed and time of the timer are adjustable in this mode.

1.4 Turn off the device only after pressing "STOP".

2. Remote control:

2.1 When plugged in, press the "  " -button and the display will show:

"Default time: 10 minutes. Speed: 00."

2.2 Press the "  /  "-button, the unit starts in P1 mode by default.

(View the program overview of "P1" for details)

2.3 If you want to enjoy P2,- P3,- P4 or P5 mode, press the "  /  "-button to exit P1 mode before selecting P2,- P3,- P4 or P5 mode.

Then press the P2,- P3,- P4 or P5 mode first, and then press the "  /  "-button to start the mode.

(Please refer to the program overview of "P1", "P2", "P3", "P4" and "P5" for details)

2.4 Turn off the unit only after pressing the "  /  " button.

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

- ⚠ Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.
- Turn off and unplug the product before cleaning.
- Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its life, do not throw the product away with the normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl
- We provide one year warranty on our products. The customer must provide proof of purchase and the warranty must be validated within a period of one year from the date of purchase. Please note the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Disclaimer

- Changes subject to change; specifications are subject to change without notice.

Troubleshooting guide

Common problems	Causes	Solutions
The massager does not work	1. Disconnect power or loose plug causing poor contact. 2. No power or power cord is damaged 3. The power button is off, or the timer has been reset. 4. The motor auto shuts-off due to the device being on too long activating overheating protection. 5. Didn't use original adapter.	1. Remove plug from outlet and re-insert into power source. 2. Try another outlet. Check the power cord, replace if damaged. 3. Turn on the power button. 4. Use the massager after the motor cools off naturally. 5. Use original adapter provided in the package only.
The massage intensity is too high or too low	The incorrect intensity level for your preference was selected.	The massager has 20 speed levels, choose the speed you prefer. We recommend low speed to begin with.
The time and the speed cannot be adjusted	In the P2,- P3,- P4- and P5-modes, the time of the timer and speed cannot be adjusted and no button is available except the Start/stop-button.	Select P1-mode.
Rotation doesn't work	The rotation mode of the massage heads is different for each mode.	Select the desired mode and activate the desired mode: P1, P2, P3, P4 or P5.

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© OnlineGoods