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Introduction

Thank you for choosing our product! Always follow the safety instructions for safe use. If there are any questions about the product please feel free to contact us via email: info@good-4you.nl

Intended Use

- This product is intended only for the purposes described in the instructions for use.
- The supplier is released from liability for any injury or damage caused by using the product beyond its intended use.
- Always observe the following precautions when using "the Foot Massager", hereinafter referred to as "the product".

Warnings and Safety Instructions

Make sure that the packaging material does not get into the hands of children. They may choke on it! Read the instructions for use carefully before use and keep them carefully!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully. Be aware of the risks and consequences associated with misuse of the product. Misuse of the product may result in personal injury and damage to the product or products connected to it.
- Use the product only for the purposes indicated in this manual.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult Customer Service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Do not leave the product unattended when in use. When not in use, turn the product off.

Persons

- The product is not suitable for people (including children) with physical or mental disabilities, or people who have no experience with or knowledge of the product. Use of the product is permitted only when a supervisor is present and responsible for the safety of the user. Use and store the product out of the reach of children. We recommend using the product only if you are 16 years of age or older.
- Do not use the product if you: have an electronic implant such as a pacemaker; have conditions on your foot(s) such as warts, swimmer's eczema, rashes, or open wounds; have a malignant tumor, heart disease, osteoporosis, or a broken spine; have a cold with a fever or varicose veins, thrombosis, jaundice, diabetes, nervous disorders, acute inflammation, skin disease, broken skin, or pain without obvious or immediate cause.
- Do not use the product with more than one person at a time. Do not use the product near your eyes or other sensitive areas. Do not use the product when you are in bed or sleeping.
- Do not use the product when you are driving a vehicle or operating machinery. Do not use the product vertically to avoid accidents or injuries to the turbine, malfunctions, etc.
- If you lend the product to someone else, make sure you give the product with the manual together.
- Use the product on a level surface. If it is not stable, the body will not be balanced and an accident can be caused.
- The product is special for the legs, feet, calves, etc., do not use the waist, neck and other parts outside the purpose of use.

- Prolonged use of the product may cause overuse. After 15 minutes of use, insert a break of at least 15 minutes before using the product again.
- Do not leave the product running unloaded or reverse the direction of rotation of the product. Friction can be harmful to the skin, so do not use it directly.
- Do not place heavy objects on the product.
- Do not allow the product to come into contact with pointed or sharp objects.

Medical Notes

- The product is not a substitute for treatments from a doctor or physical therapist.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment, or are using medical equipment.
- To avoid excessive over stimulation of the muscles and nerves, it is recommended that you do not use the product for more than 15 minutes at a time.
- If you experience pain or discomfort while using the product or find the massage unpleasant, discontinue use and consult your physician.

Electra

- Do not disassemble, repair or restore the product. Doing so may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to prevent damage and personal injury.
- Do not use the product near a bathtub, shower, swimming pool, or near a tub or laundry tub filled with water.
- Do not allow the product to come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service department.
- Do not touch the product with wet hands.
- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This may cause a fire hazard, electric shock, personal injury or damage to the product.

Other

- Be sure to press the buttons briefly to avoid confusing the product.
- Do not insert objects into or through the openings of the product. Always ensure that the moving parts have room to move and that nothing is jammed.

Receiving Shipment

- Check that the contents of the package match the packing slip when you receive the shipment. Notify your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and give the supplier a detailed description of the damage. Keep the damaged packaging to substantiate your complaint.

Package Contents

- Good4You Shiatsu Foot Massager - Model XB-8584

Attached accessories

- 1x Remote Control
- 1x CR2025 Battery
- 1x Power Adapter
- 1x User Manual

Instructions for use

- ⚠ Prolonged use of the product may lead to overheating. After 15 minutes of use, take a break of at least 15 minutes before using the product again and do not use it more than 3 times a day.
- ⚠ Do not cover the product when it is on. Never place pillows or blankets on or against the unit. This may cause a fire hazard, electric shock, personal injury or damage to the product.
- ⚠ Keep the product away from moisture, heat, areas where the temperature is relatively high, and heat sources such as ovens.

Health benefits

- Foot reflexology has become increasingly popular in health care as people place more emphasis on wellness. Daily activities such as long walks or running errands can be tiring, and massaging your feet at home can provide a sense of relief.
- According to traditional Chinese medicine, internal organs and legs are closely connected. The human ankle has more than 60 reflex points that, when stimulated, can improve blood circulation and organ function. Massaging your feet before going to bed reduces fatigue, calms your nerves and improves sleep, contributing to a better overall quality of life.


Application of the Foot Massager

- Foot reflexology reflexes are responses to external stimuli from our bodies. When the tissues and organs of our body are not functioning properly in the corresponding areas on our feet, air bubbles, small nodules or a grainy sensation may occur, and these areas may be sensitive. This pain is transmitted through the afferent nerve to the central nervous system, and this discomfort can affect various parts of the body.
- Treatment for headaches, abdominal pain, migraines, neck and shoulder pain, fatigue and blood pressure problems.

Product Theory of Foot Massage Device

- Reflexology is based on traditional medicine, meridian theory, magnetic sciences and health legislation concerning reflexology. It combines modern electronic technology with medical research results.

Principles of neural reflexes

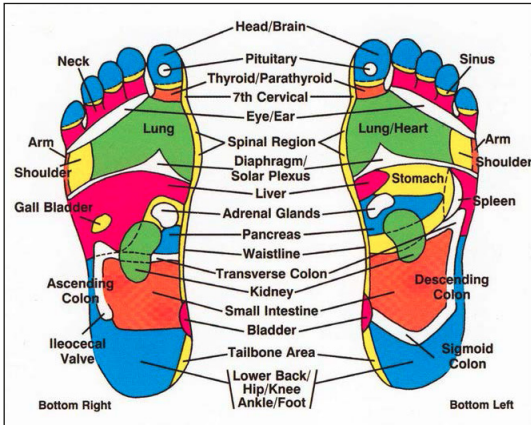- The reflex is a human response to external stimuli. Abnormal body tissues and organs can cause problems at the reflex points, such as air bubbles, grains of sand, rod-like locks and small nodules. Pain is sent by the afferent nerves to the central nervous system, which then sends a series of nerve impulses to the body tissues and organs. This leads to hormonal regulation and increases the body's ability to adapt to imbalances.

Features of the Good4You® Foot Massager.

- Complete massage with focus on key foot therapy points:** The device provides comprehensive treatment by targeting specific pressure points in the feet, increasing the effectiveness of therapy.
- Simultaneous dual massage on both feet:** Research has shown that simultaneous stimulation of pressure points provides a benefit that is 3 to 5 times greater than the traditional method.
- Acupressure points:** The foot massager works on each area of the foot to treat different types of pathologies according to the zones defined in foot reflexology.
- Ease of use at any time:** It's like bringing home a professional, with the advantage that you can perform massage sessions after a shower, before bed or while watching television. Calibri
- Unique design and shape:** For perfect distribution and adaptation to any foot size and pressure points. The control panel is easy to access and operate, allowing you to take advantage of all the functions of your foot massager.
- Advanced microcontroller:** Ensures optimal performance during use while reducing power consumption.
- Customizable:** Can be adjusted to your personal preferences at any time during your massage session.

Frequently Asked Questions

1. **When should I not use the foot massager?**

It is not recommended to use the device 30 minutes before and 1 hour after eating. This may reduce blood circulation in the gastrointestinal area and affect digestion.

2. **How long is recommended to use the foot massager?**

It is recommended to use the device for 30 minutes every session.

3. **What care should I have after each massage?**

It is recommended to drink a cup of warm water to promote blood flow and circulation. Tea, wine and other alcoholic beverages should be avoided.

4. **Is the use recommended for women and men?**

Yes, use is recommended for both women and men. It can contribute to a significant improvement in your health. By properly stimulating the pressure points in the feet, the adrenal glands can be stimulated to release adequate amounts of hormones into the body, which vitalizes the skin, reduces pigmentation and speeds up metabolism to reduce body weight. In women, it can help normalize menstrual periods.

What reactions and results may occur after using the foot massager?

After using the device 5-10 times, you may notice the following improvements in your overall health:

1. **Improved sleep quality:** During massage, you feel sleepy and may want to sleep. This reaction indicates that the physiological functions of the body are being adjusted.
2. **Increased perspiration:** This helps to remove toxins from the body and improves blood circulation.
3. **Prevent constipation:** It promotes proper functioning of the digestive system by normalizing bowel movements.
4. **Patients with back or shoulder pain may experience some discomfort during foot massage:** This indicates that damage to the body is being repaired and the pain will disappear completely over time.
5. **Patients with lymphatic problems may experience inflammation in the feet:** Keep performing the massage so that the fluids can circulate, and the inflammation disappears naturally.
6. **Feeling warm:** You may experience a feeling of warmth due to the increased blood circulation from the stimulation of the reflex zones in the feet.

Common diseases and foot therapy to treat them

Headache and migraines

Massage the indicated area to reduce the headache, you can also do it preventively to maintain a healthy state.

Shoulder pain

Massage the indicated area to relieve pain and pressure on the shoulders to promote blood circulation in the upper body.

Cervical pain

Neck and shoulder pain are usually related, so it is recommended to simultaneously massage the areas indicated in both images. It is recommended to massage the designated area in the morning.

Stomach pain and upset

Massage the indicated area to relieve stomach pain and discomfort such as gastritis, vomiting, etc.

Fatigue

Massage the targeted area to adjust neural reflexes, adjust endocrines, improve circulation and relieve muscle tension.

Eliminate constipation

Massage the indicated area to stimulate a bowel movement, protect the liver, kidneys, and other important organs.

Weight reduction

By massaging the indicated area, it will directly influence the digestive system, improving its functioning with benefits in weight reduction.

Detoxification

By massaging the indicated areas, the elimination of toxins through urine is promoted to improve kidney function. The inner part of the heel is related to the bladder. The lateral part of the heel is at the top of the lymph nodes. The big toe line (Hallux) is related to the lower lymph nodes.

Lack of sleep, insomnia

Massage the indicated areas, this will allow you to fall asleep and reduce stress.

Colds

The indicated areas are related to: nose, lungs, bronchi, larynx and other respiratory reflex zones. Massaging this area will help relieve symptoms, improve circulation, and achieve a physical state of balance.

Stress and pressure

Everyday stress can be seen reflected in the big toe with a feeling of hardness to the touch. By massaging the finger, it will be soothed again and the pituitary gland will be stimulated, which will release hormones regulating anxiety levels, controlling depression, and providing mental relaxation.

Dysmenorrhea

Massage the indicated areas until you feel comfortable and the area is soothed. It is recommended to massage the indicated point to treat dysmenorrhea from the first day of menstruation so that the process runs normally and without pain.

Climacteric syndrome (Menopause)

Massaging the indicated points will allow you to alleviate symptoms related to menopause.

Infertility and menstrual problems

The massage of the indicated areas improves the function of the ovaries, the meridian located in the lower part of the heel is the point related to the uterus, a massage at this point can relieve pain.

Start-up instructions


1. Plug the power cord into the socket.
2. Press the "⏻" button to turn on the device. The corresponding LED-indicator will turn red.
3. Press the "Auto" button to switch the massage mode between Programs 1, 2, 3, and 4:
 - a) Program 1: General massage
 - b) Program 2: Massage focused on head, neck, spine, and shoulders
 - c) Program 3: Massage focused on stomach, intestines, liver, and other organs
 - d) Program 4: Massage focused on legs, urinary system, and reproductive system
 The corresponding LED indicator will turn red.
4. Press the "Manu" button for a personalized massage; the massage mode alternates between movements for the toes, arches, and soles. Use the "Speed" button to adjust the speed of the massage. The corresponding LED indicator will turn red.
5. Press the "⏻" button to select the manual mode and adjust the speed and direction as desired. The corresponding LED-indicator will turn red.
6. Press the "⏻" button and hold it for 2 seconds to turn off the device. Unplug the power cord from the socket after use.

- Auto-shutdown: If no massage mode is selected within 5 minutes, the device will automatically turn off.
- During the massage: Press the "⏻" button to pause the operation and press again to resume.

⚠ During operation: never pull the power cord to turn off the unit.

- For people with sensitive feet, the device may feel a bit rough at first. It is recommended to lift your feet slightly at the beginning. Over time, your feet will adapt to the device after several uses. We are confident you will eventually love it.

Operation panel

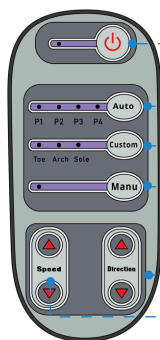On/Off:
Starts or ends the massage.
After 15 minutes, it stops automatically.


Automatic stands:
4 Massage Programs

Personalized Massage
3 Massage zones: Toes - Arches - Soles

Manually Adjust Speed and Direction

Adjust Direction

Adjust Speed

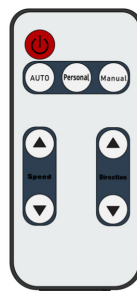All associated functions can be set directly with the remote control.

Technical data

Dimensions	56 x 33 x 25 cm
Weight	4,9 kg
Power	40W
Power - Frequency - Output	120V - 50/60Hz - 1.5A
Timer	15 min

Maintenance and Cleaning

Maintenance

- Turn off the product when not in use and unplug it.
- Keep the product dust-free.
- Make sure the product is not in contact with sharp objects that may damage the product.
- Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

Cleaning

- ⚠ Do not immerse the product in water or other liquids and keep it away from solvents and harsh cleaning agents.
- Turn off and unplug the product before cleaning.
 - Clean the base and control panel of the product with a damp cloth dipped briefly in detergent and wipe the product dry with a dry, soft cloth. Do not use abrasives or thinners to clean the product.

Disposal and Recycling

- Think about the environment and contribute to a cleaner living environment! At the end of its service life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Service and Warranty

- If you need service or information regarding your product, please do not hesitate to contact us via email: info@good-4you.nl. You will receive a response within 24 hours.
- We provide two years warranty on our products. The customer must be able to show proof of purchase and the warranty must be made valid within a period of two years from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects known to the customer at the time of purchase; for defects caused by the customer's own fault; for damage caused by third parties.

Extend warranty for FREE

1. Scan the QR code below to find Good4You's Facebook page and then 'Like' it

2. Send a message with "Warranty" to us via Messenger or email info@good-4you.nl to extend the warranty from 2 years to 4 years.

Troubleshooting

Problem	Solution
A clattering sound is heard in the movement.	This is perfectly normal and is caused by the operation of the product. It is a normal status.
The device suddenly turns off	The device automatically turns off after 15 minutes. Continue use by pressing the power button. If the device has been used for 30 minutes, turn the device back on after resting for 15 minutes.
	Overheating protection is activated. The appliance automatically turns off for safety. Let the appliance rest for 15- 30 minutes before using it again.
The device will not start	Check if the device has been used for more than 15 minutes. If so, turn the unit back on after resting for 15 minutes. In warmer environments, it may take longer for the device to cool down.
	Check if the adapter is plugged in.

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